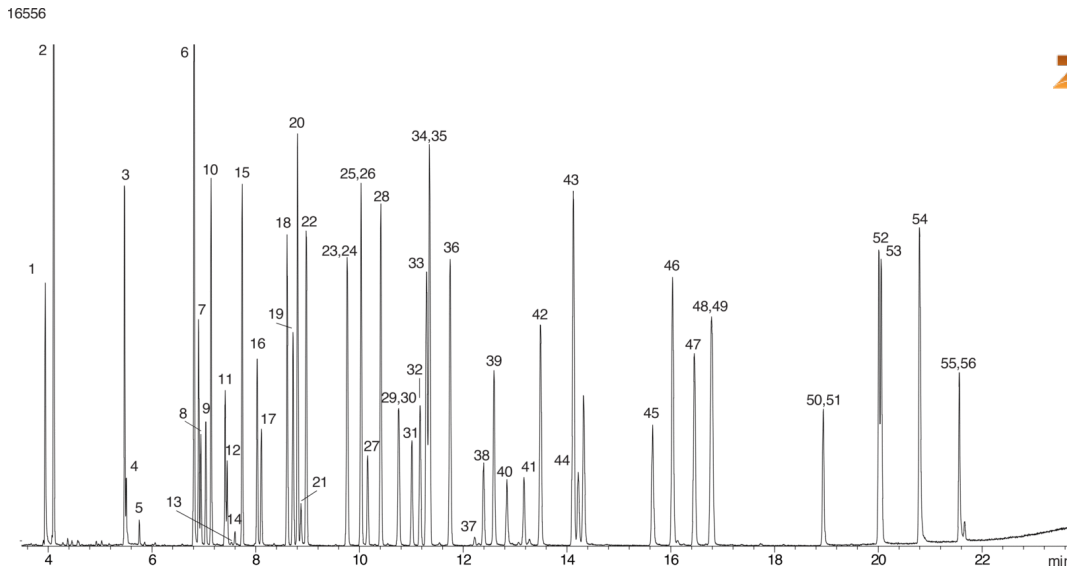


Pesticides by MSD on ZB-5MSi

Column: Zebron™ ZB-5MSi, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea
Phase: 5% Phenyl 95% Dimethylpolysiloxane
Dimensions: 30 meters x 0.25 mm x 0.25 µm
Order No: 7HG-G018-11
Oven Profile:

Carrier Gas: Constant Flow Helium, 1.3 mL/min
Injection: Split 20.000000000:1 0 µL @ 250°C
Detection: Refractive Index (180°C)
Analyst Note: Anaytes are 10ppm in Acetone

**Products used in this application:**

Pesticides by MSD on ZB-5MSi**ANALYTES:**

- | | |
|-------------------------------|-----------------------------|
| 1 Methamidophos | 41 Methidathion |
| 2 Dichlorvos | 42 Tetrachlorvinphos |
| 3 Mevinphos | 43 Tokuthion |
| 4 Acephate | 44 Profenofos |
| 5 Trichlorfon | 45 Fensulfothion |
| 6 Omethoate | 46 Ethion |
| 7 Thionazin | 47 Bolstar |
| 8 Demeton | 48 Carbophenothion |
| 9 Demeton-S-methyl | 49 Famphur |
| 10 Ethoprop | 50 Phosmet |
| 11 Naled | 51 EPN |
| 12 Dicrotophos | 52 Azinphos methyl |
| 13 Monocrotophos | 53 Leptophos |
| 14 Sulfotep | 54 Azinphos-ethyl |
| 15 Phorate | 55 Pyrazophos |
| 16 Demeton | 56 Coumaphos |
| 17 Dimethoate | 57 Dioxathion |
| 18 Terbufos | |
| 19 Fonofos | |
| 20 Diazinon | |
| 21 Phosphamidon | |
| 22 Disulfoton | |
| 23 Dichlofenthion | |
| 24 Phosphamidon | |
| 25 Methyl parathion | |
| 26 Chlorpyrifos-methyl | |
| 27 Tolclofos-methyl | |
| 28 Fenchlorphos | |
| 29 Fenitrothion | |
| 30 Pirimiphos-methyl | |
| 31 Malathion | |
| 32 Aspon | |
| 33 Fenthion | |
| 34 Chlorpyrifos | |
| 35 Parathion | |
| 36 Trichloronate | |
| 37 Chlorfenvinphos | |
| 38 Merphos | |
| 39 Chlorfenvinphos | |
| 40 Crotoxyphos | |

