


**Endocrine disruptors by MSD on ZB-5MSi****Column:** Zebtron<sup>™</sup> ZB-5MSi, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G018-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1.2 mL/min**Injection:** Split 40.0000000000:1 0 µL @ 250°C**Detection:** Refractive Index (180°C)**Analyst Note:** Anaytes are 50ppm in Acetone**Products used in this application:**

## Endocrine disruptors by MSD on ZB-5MSi

### ANALYTES:

- 1 Dimethyl phthalate
- 2 Diethyl phthalate
- 3 Atraton
- 4 Simazine
- 5 Prometon
- 6 Atrazine
- 7 Propazine
- 8 Dipropyl phthalate
- 9 Terbutylazine
- 10 Secbumetone
- 11 Simetryn
- 12 Ametryn
- 13 Prometryn
- 14 Terbutryn
- 15 Dibutyl phthalate
- 16 4,4'-DDE
- 17 4,4'-DDD
- 18 Di-n-hexyl phthalate
- 19 4,4'-DDT
- 20 Dicyclohexyl phthalate
- 21 bis(2-Ethylhexyl)phthalate
- 22 Di-n-octyl phthalate
- 23 Ethinyl estradiol

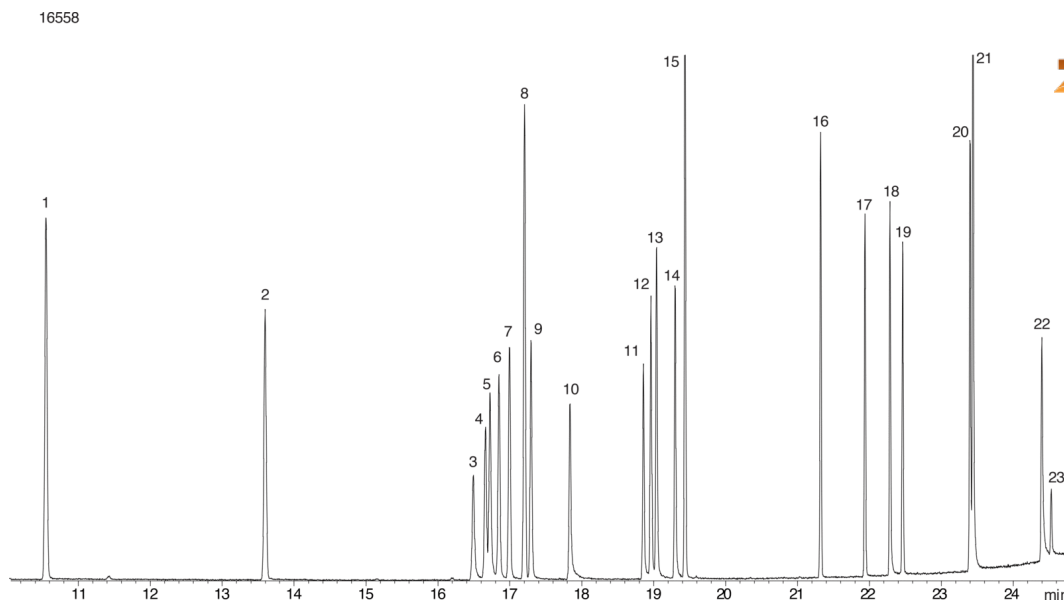