


Benzodiazepines by MSD on ZB-5MSi**Column:** Zebron[™] ZB-5MSi, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G018-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1.5 mL/min**Injection:** Splitless 0.000000000:1 0 μ L @ 250°C**Detection:** Refractive Index (180°C)**Analyst Note:** Anaytes are 5ppm in methanol**Products used in this application:**

Benzodiazepines by MSD on ZB-5MSi

ANALYTES:

- 1 Oxazepam
- 2 Lorazepam
- 3 Diazepam
- 4 Nordazepam
- 5 Midazolam
- 6 Temazepam
- 7 Flunitrazepam
- 8 Bromazepam
- 9 7-Aminoflunitrazepam
- 10 Prazepam
- 11 Lormetazepam
- 12 Flurazepam
- 13 Estazolam
- 14 Alprazolam
- 15 Triazolam

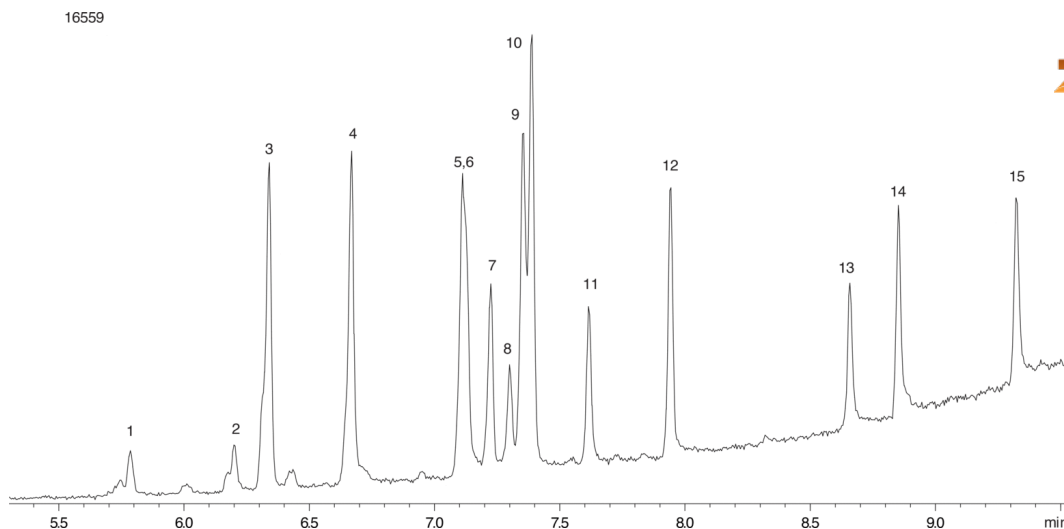