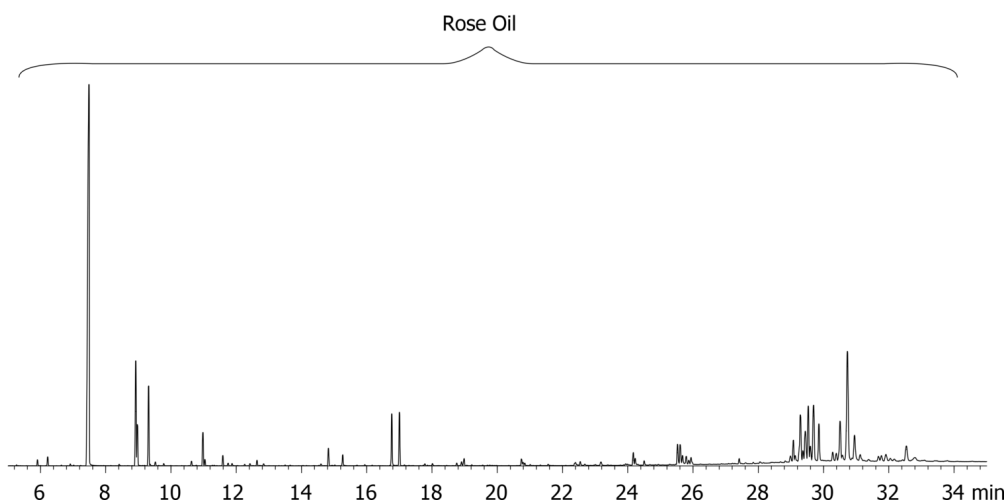


Rose Oil by MSD (Straight liner) on ZB-XLB**Column:** ZebtronTM ZB-XLB, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** Proprietary XLB Phase**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G019-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1.2 mL/min**Injection:** Split 25.0000000000:1 0 μ L @ 250°C**Detection:** Refractive Index (230°C)**Analyst Note:** Sample is 1% Rose oil in ethyl acetate

16655

**Products used in this application:****ANALYTES:****1** Rose oil