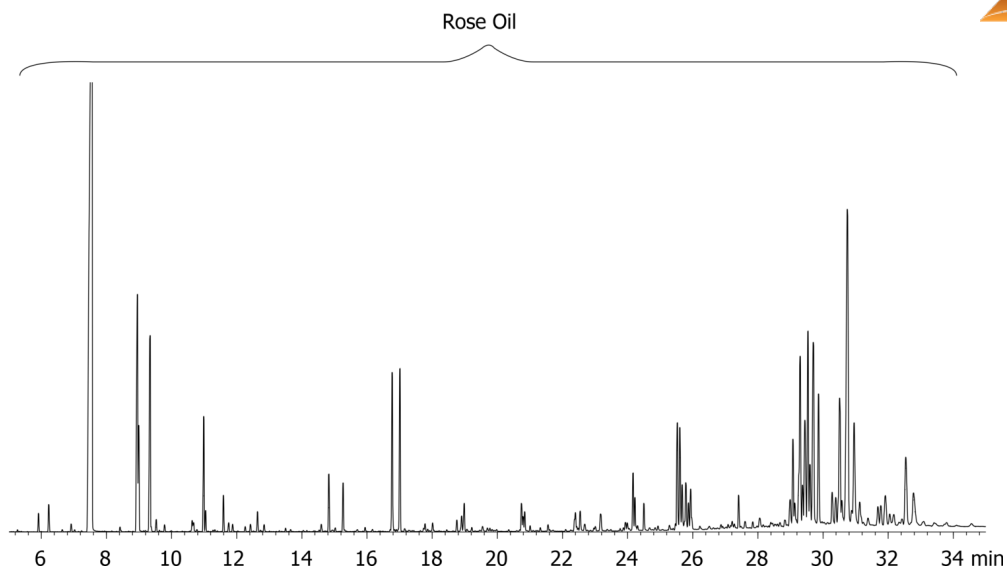


**Rose Oil by MSD (Gooseneck liner) on ZB-XLB****Column:** Zebron<sup>™</sup> ZB-XLB, GC Cap. Column 30 m x 0.25 mm x 0.25  $\mu$ m, Ea**Phase:** Proprietary XLB Phase**Dimensions:** 30 meters x 0.25 mm x 0.25  $\mu$ m**Order No:** 7HG-G019-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1.2 mL/min**Injection:** Split 25.0000000000:1 0  $\mu$ L @ 250°C**Detection:** Refractive Index (230°C)**Analyst Note:** Sample is 1% Rose oil in ethyl acetate

16656

**Products used in this application:****ANALYTES:****1** Rose oil