


**ASTM D5501: Determination of Denatured Bioethanol on ZB-Bioethanol****Column:** Zebron™ ZB-Bioethanol, GC Cap. Column 15 m x 0.25 mm x 1.00 µm, Ea**Phase:** Proprietary Bioethanol Phase**Dimensions:** 15 meters x 0.25 mm x 1 µm**Order No:** 7EG-G020-22**Oven Profile:****Carrier Gas:** Constant Flow Helium, 25 cm/sec**Injection:** Split 50.0000000000:1 0 µL @ 300°C**Detection:** Refractive Index (300°C)**Analyst Note:** Instrumentation: Shimadzu GC-2010**Products used in this application:**

16667

**ANALYTES:**

- 1 Methanol
- 2 Ethanol

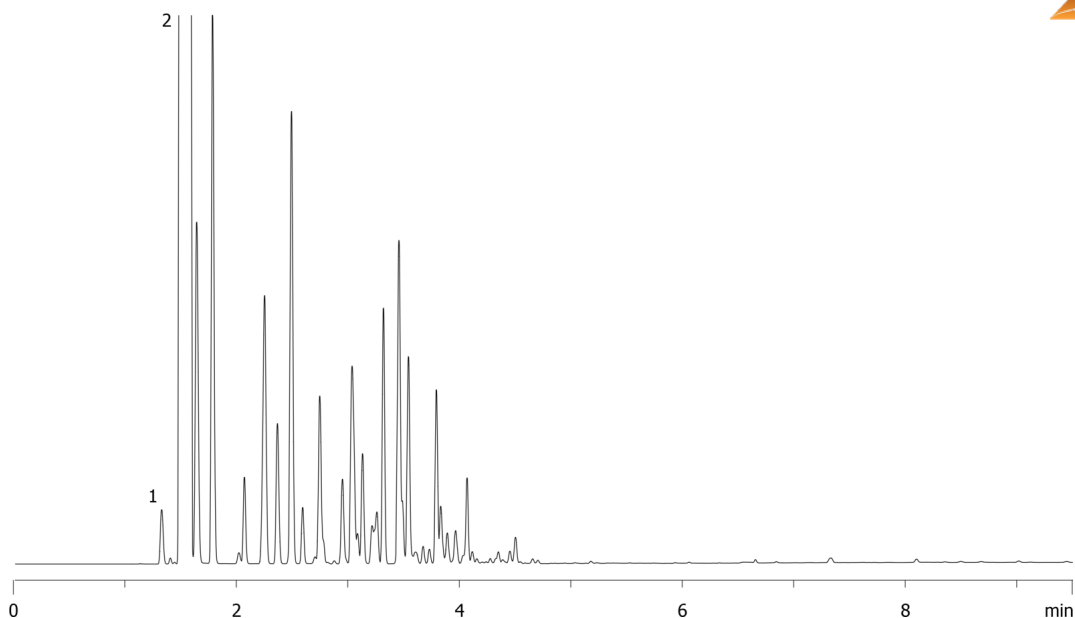