

# HPLC Application

ID No.: 16677

## Mixtures of Acids, Neutrals and Bases in 10mM Ammonium Acetate (pH 4.8) on Gemini NX-C18 5µm

**Column:** Gemini® 5 µm NX-C18 110 Å, LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-4454-E0

**Elution Type:** Gradient

**Eluent A:** 10mM Ammonium Acetate pH 4.8

**Eluent B:** ACN

| Gradient | Step No. | Time (min) | Pct A | Pct B |
|----------|----------|------------|-------|-------|
| Profile: | 1        | 0          | 95    | 5     |
|          | 2        | 8          | 20    | 80    |
|          | 3        | 0          | 20    | 80    |

**Flow Rate:** 1.5 mL/min

**Col. Temp.:** 25 °C

**Detection:** Refractive Index @ 0.0000000000 (ambient)

**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Gemini-NX C18 4 x 3.0mm, 10/Pk Part No.: AJ0-8368


- Holder Part No.: KJ0-4282



Products used in this application:



16677

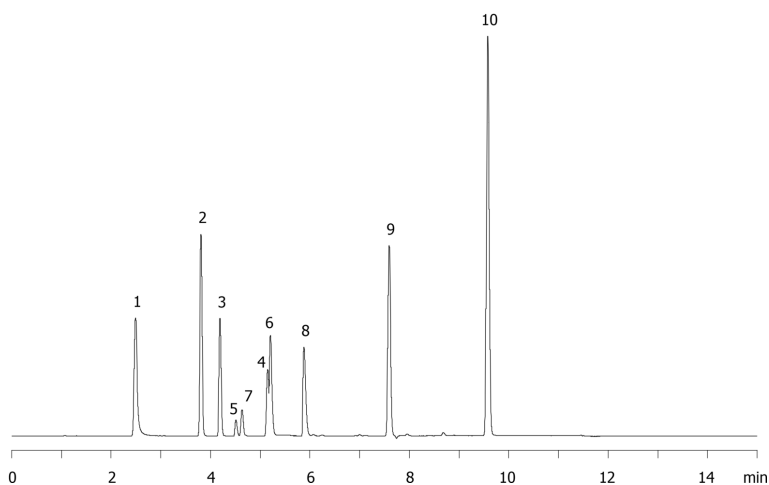

## Mixtures of Acids, Neutrals and Bases in 10mM Ammonium Acetate (pH 4.8) on Gemini NX-C18 5µm

### ANALYTES:

- 1 Pyridine
- 2 Quinidine
- 3 Sulfathiazole
- 4 Triprolidine
- 5 Benzyl alcohol
- 6 Nortriptyline
- 7 Phenol
- 8 3-Methyl-4-nitrobenzoic acid
- 9 Methyl salicylaldehyde
- 10 Hexanophenone

