

## Methoxyphenamine on Lux 3 um Cellulose-1

**Column:** Lux® 3 µm Cellulose-1, LC Column 250 x 4.6 mm, Ea

**Dimensions:** 250 x 4.6 mm ID

**Order No:** 00G-4458-E0

**Elution Type:** Isocratic

**Eluent A:** 0.1% DEA in Hexane

**Eluent B:** 0.1% DEA in Isopropanol

| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B |
|-------------------|----------|------------|-------|-------|
| 1                 | 0        | 80         | 20    |       |

**Flow Rate:** 1 mL/min

**Col. Temp.:** 25 °C

**Detection:** Refractive Index @ 0.0000000000 (ambient)

**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.


- SecurityGuard Cartridges, Lux Cellulose-1, 4 x 3.0mm, 10/Pk Part No.: AJ0-8403

- Holder Part No.: KJ0-4282

16719



Products used in this application:

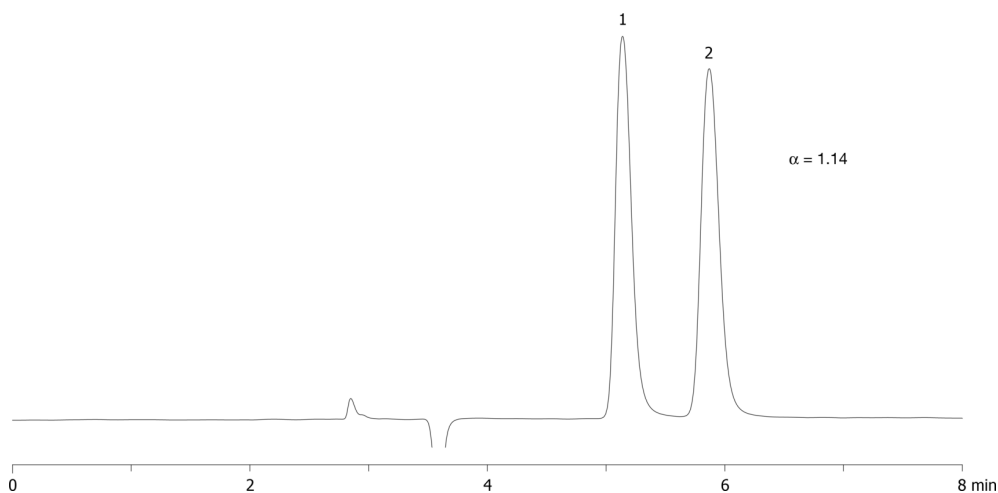

### ANALYTES:

- 1 Methoxyphenamine
- 2 Methoxyphenamine

