


Methadone by MSD on a ZB-5MS**Column:** Zebtron[™] ZB-5MS, GC Cap. Column 10 m x 0.18 mm x 0.18 µm, Ea**Phase:** 5% Phenyl-Arylene 95% Dimethylpolysiloxane**Dimensions:** 10 meters x 0.18 mm x 0.18 µm**Order No:** 7CD-G010-08**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1.2 mL/min**Injection:** Split 10.0000000000:1 0 µL @ 250°C**Detection:** Refractive Index (320°C)**Analyst Note:** Analytes are 24ppm
Sim ions:
(Quant ion*) D9-Methadone (*78,303): Methadone (*72,294,223)**Products used in this application:****ANALYTES:**

- 1 Methadone-d9
- 2 Methadone

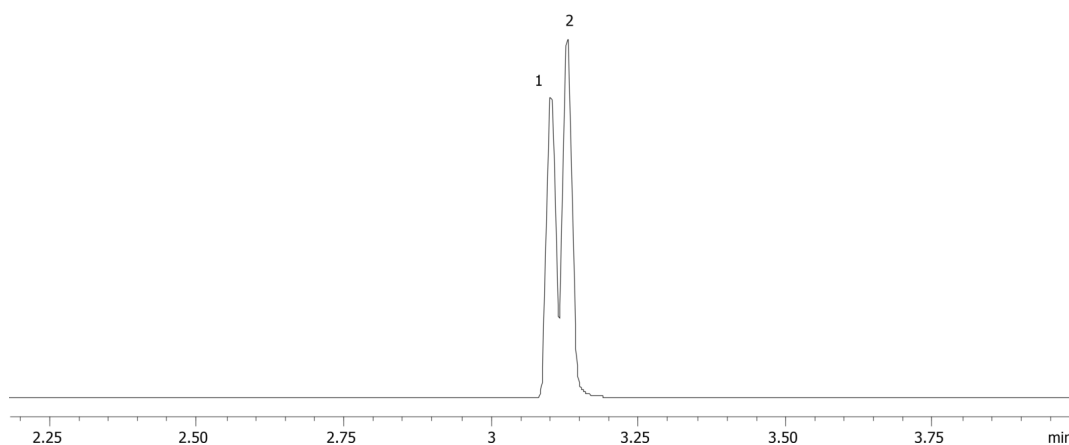