

Sterols from Olive Oil by FID on ZB-5HT**Column:** Zebtron™ ZB-5HT, GC Cap. Column 30 m x 0.25 mm x 0.10 µm, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.1 µm**Order No:** 7HG-G015-02**Oven Profile:****Carrier Gas:** Constant Flow Helium, 2 mL/min**Injection:** Splitless 0.000000000:1 0 µL @ 350°C**Detection:** Refractive Index (350°C)**Analyst Note:** Sample Prep: Saponification, Solid Phase Extraction (SPE), and Derivatization.

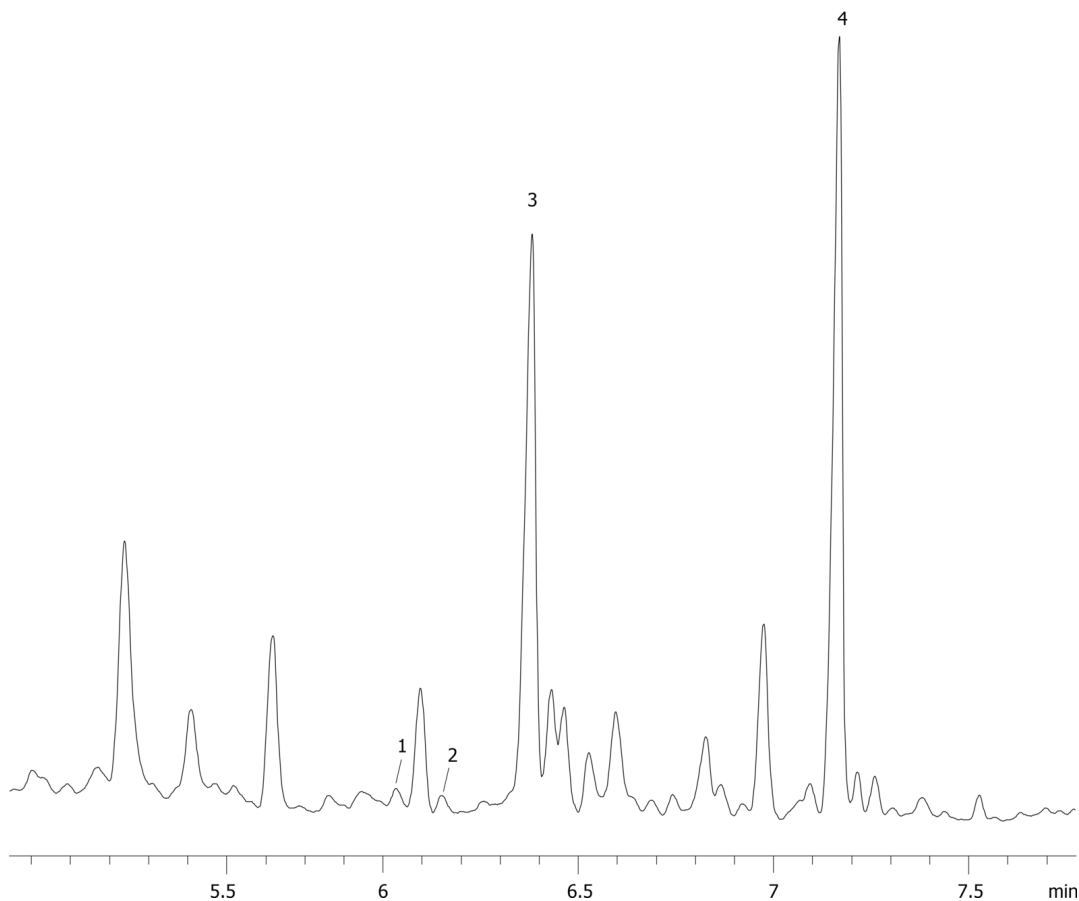
Recommended Accessories:

Recommended Liner: Zebtron PLUS Single Taper Z-Liner, 4 mm ID

Liner Part No.: AG2-0A13-05 (for Agilent systems)

**Products used in this application:**

16753

**ANALYTES:**

- 1 Campesterol
- 2 Stigmasterol
- 3 b-Sitosterol
- 4 Betulin

