


Sterols by FID on ZB-MR 2**Column:** Zebron™ ZB-MultiResidue™-2, GC Cap. Column 30 m x 0.25 mm x 0.20 µm, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.25 mm x 0.2 µm**Order No:** 7HG-G017-10**Oven Profile:****Carrier Gas:** Constant Flow Helium, 2 mL/min**Injection:** Splitless 0.000000000:1 0 µL @ 300°C**Detection:** Refractive Index (300°C)**Analyst Note:** Sterols derivatized via BSTFA:TMCS**Products used in this application:**

16754

**ANALYTES:**

- 1 Cholesterol
- 2 Campesterol
- 3 Stigmasterol
- 4 b-Sitosterol
- 5 Betulin

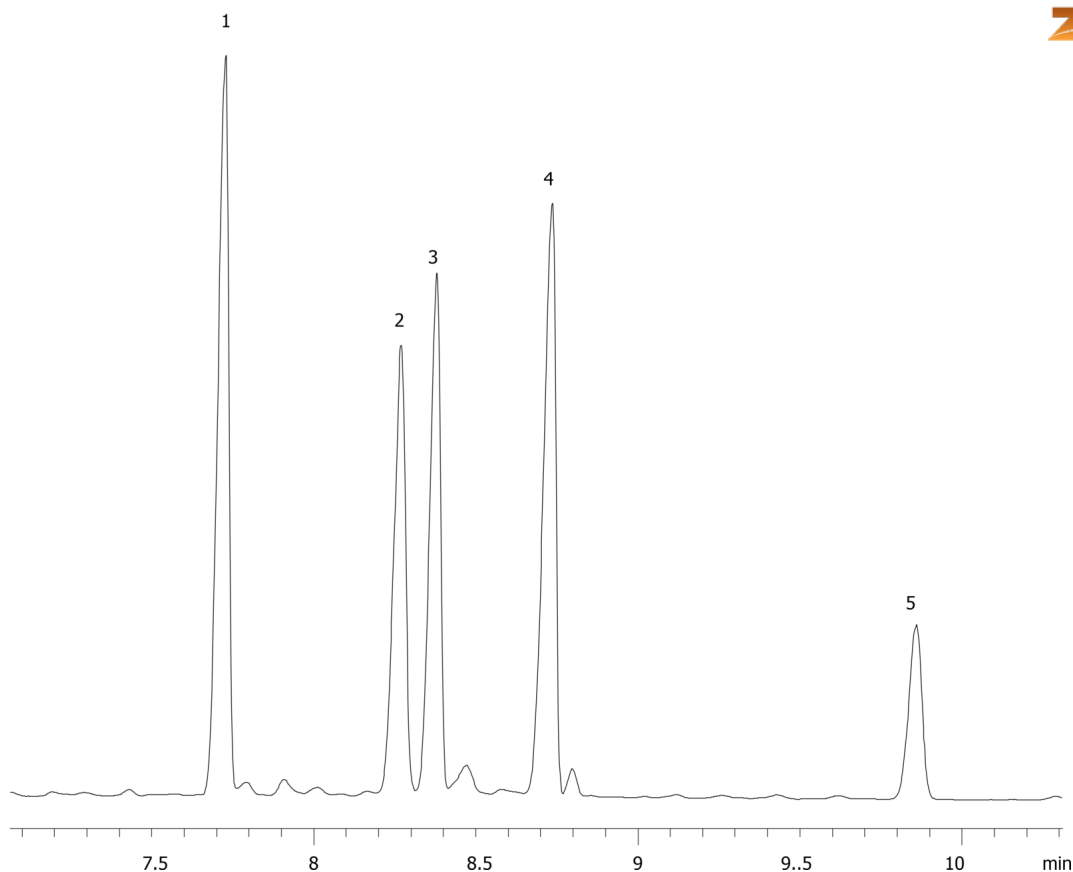