


EPA Method 8081A on ZB-35 (Carrier Gas: Helium)**Column:** Zebron[™] ZB-35, GC Cap. Column 30 m x 0.32 mm x 0.25 µm, Ea**Phase:** 35% Phenyl 65% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.32 mm x 0.25 µm**Order No:** 7HM-G003-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 2.6 mL/min**Injection:** Splitless 0.000000000:1 0 µL @ 250°C**Detection:** Refractive Index (335°C)**Analyst Note:** Analytes are 50 ppb in isooctane**Products used in this application:**

EPA Method 8081A on ZB-35 (Carrier Gas: Helium)

ANALYTES:

- 1 TCMX
- 2 alpha-BHC
- 3 gama-BHC
- 4 b-BHC
- 5 Heptachlor
- 6 delta-BHC
- 7 Aldrin
- 8 Heptachlor epoxide
- 9 gamma-Chlordane (trans)
- 10 alpha-Chlordane
- 11 Endosulfan I
- 12 4,4'-DDE
- 13 Dieldrin
- 14 Endrin
- 15 4,4'-DDD
- 16 Endosulfan II
- 17 4,4'-DDT
- 18 Endrin aldehyde
- 19 Endosulfan sulfate
- 20 Methoxychlor
- 21 Endrin ketone
- 22 Decachlorobiphenyl

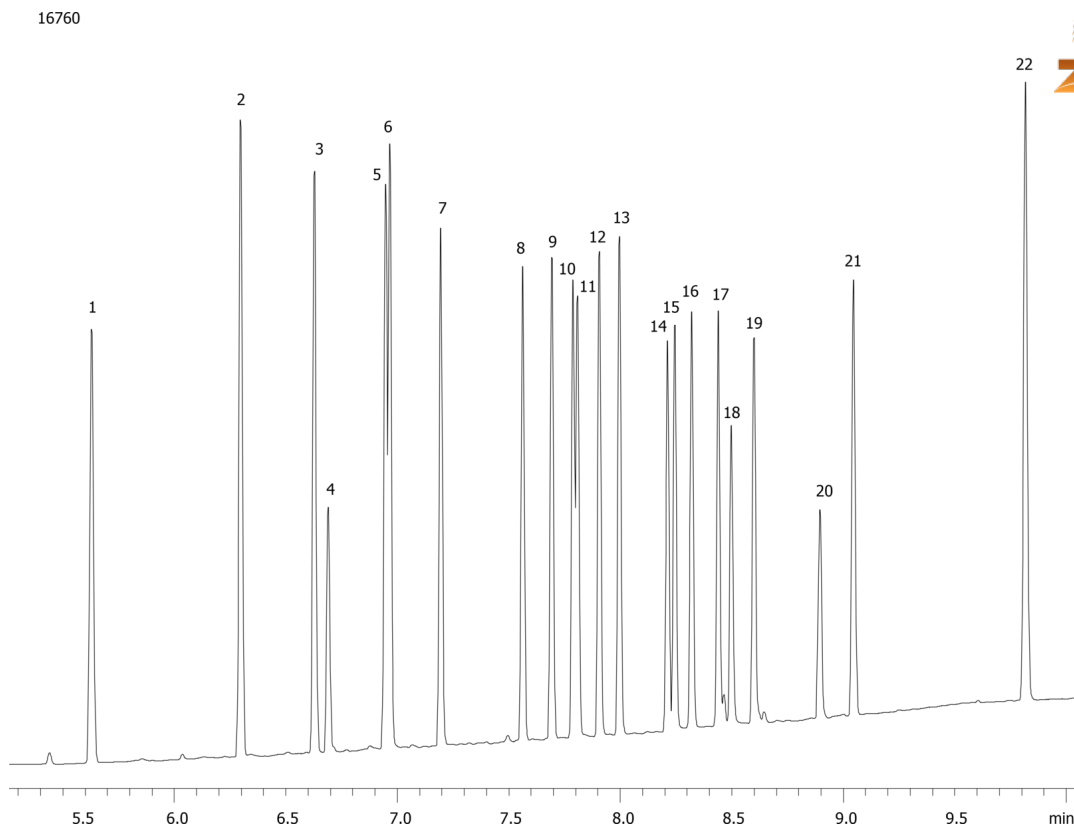