


EPA Method 552.2 on ZB-35 (Carrier Gas: Helium)**Column:** Zebron[™] ZB-35, GC Cap. Column 30 m x 0.32 mm x 0.25 μ m, Ea**Phase:** 35% Phenyl 65% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.32 mm x 0.25 μ m**Order No:** 7HM-G003-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 4 mL/min**Injection:** Splitless 0.000000000:1.0 μ L @ 250°C**Detection:** Refractive Index (300°C)**Analyst Note:** Analyte mix concentration, 2 ppm**Products used in this application:**

EPA Method 552.2 on ZB-35 (Carrier Gas: Helium)

ANALYTES:

- 1 Chloroacetic acid
- 2 Bromoacetic acid
- 3 Dichloroacetic acid
- 4 Dalapon
- 5 Trichloroacetic acid
- 6 1,2,3-Triuchloropropane
- 7 Bromochloroacetic acid
- 8 Bromodichloroacetic acid
- 9 Dibromoacetic acid
- 10 2,3-Dibromopropionic acid
- 11 Chlorodibromoacetic acid
- 12 Tribromoacetic acid

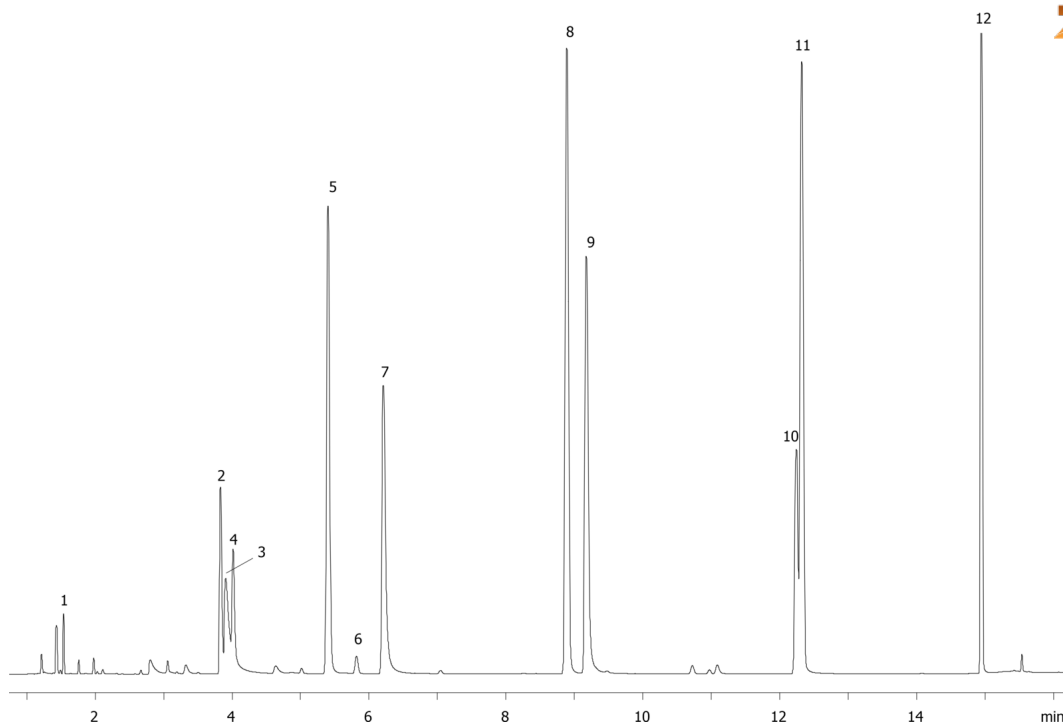