


EPA Method 8151A on MR1 (Carrier Gas: Helium)**Column:** Zebron[™] ZB-MultiResidue[™]-1, GC Cap. Column 30 m x 0.32 mm x 0.50 μ m, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.32 mm x 0.5 μ m**Order No:** 7HM-G016-17**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1.5 mL/min**Injection:** Splitless 0.000000000:1 0 μ L @ 250°C**Detection:** Refractive Index (335°C)**Analyst Note:** Analytes are 0.5 ppm in Hexane.**Products used in this application:**

16764

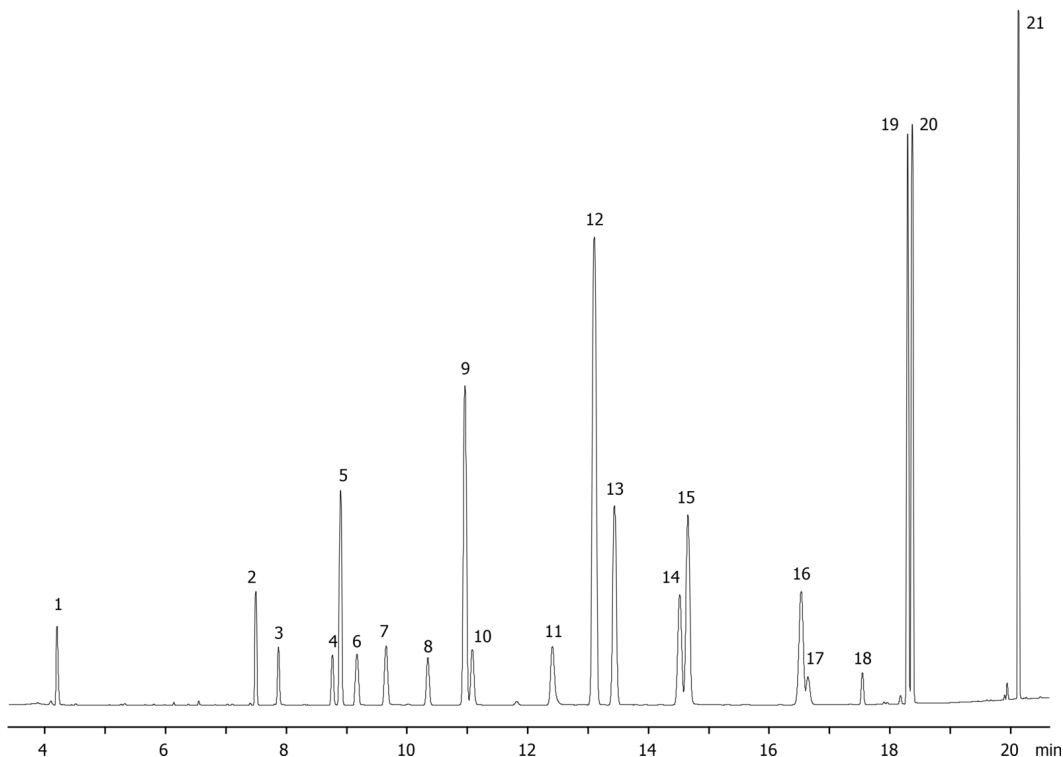

EPA Method 8151A on MR1 (Carrier Gas: Helium)

ANALYTES:

- 1 Dalapon
- 2 3,5-Dichlorobenzoic acid
- 3 4-Nitrophenol
- 4 DCAA
- 5 Dicamba
- 6 MCPP
- 7 MCPA
- 8 Dichloroprop
- 9 DBOB
- 10 2,4-D
- 11 contaminant
- 12 Pentachlorophenol
- 13 Silvex
- 14 Chloramben
- 15 2,4,5-T
- 16 Dinoseb
- 17 2,4-DB
- 18 Bentazon
- 19 Picloram
- 20 Dacthal
- 21 Acifluorfen

