


Flow Gradient Separation of a Cold Mix on Gemini C18

Column: Gemini® 5 µm C18 110 Å, LC Column 150 x 4.6 mm, Ea**Dimensions:** 150 x 4.6 mm ID**Order No:** 00F-4435-E0**Elution Type:** Gradient**Eluent A:** 10 mM Phosphate Buffer, pH 5.0**Eluent B:** Methanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	60	40
	2	3.9	60	40
	3	3.91	60	40
	4	10	60	40

Flow Rate: 0.8 mL/min**Col. Temp.:** ambient**Detection:** Refractive Index @ 0.0000000000 (ambient)**Analyst Note:** Flow gradient separation of four analytes in cold remedy on Gemini C18: Phenylephrine, Acetaminophen, Chlorpheniramine and Dextromethorphan and impurity of Maleic acid.
16773

Products used in this application:

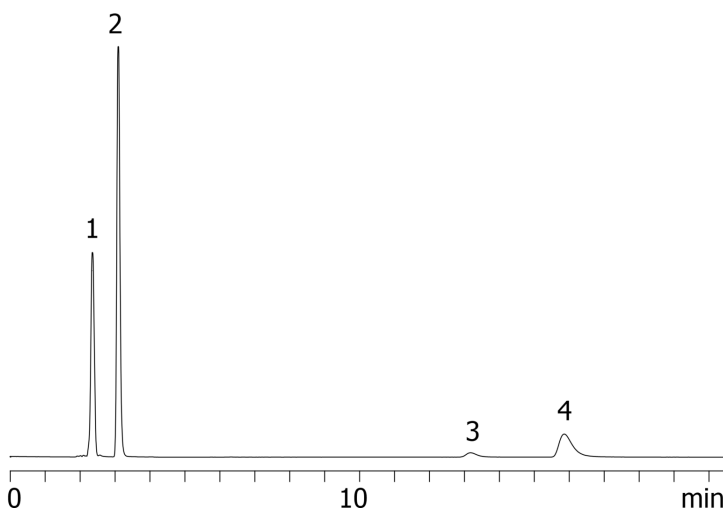

ANALYTES:

- 1 Phenylephrine
Retention Time: 2.352 min
- 2 Acetaminophen
Retention Time: 3.09 min
- 3 Chlorpheniramine
Retention Time: 13.175 min
- 4 Dextromethorphan
Retention Time: 15.854 min

