

Warfarin on Lux Cellulose-2**Column:** Lux® 5 µm Cellulose-2, LC Column 250 x 4.6 mm, Ea**Dimensions:** 250 x 4.6 mm ID**Order No:** 00G-4457-E0**Elution Type:** Isocratic**Eluent A:** 0.1% Formic acid in Hexane**Eluent B:** 0.1% Formic acid in Isopropanol

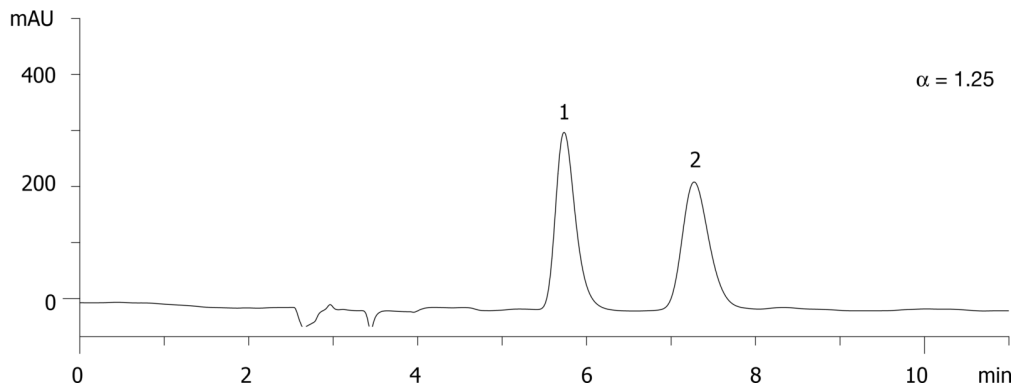
Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	20	60	40

Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** Refractive Index @ 0.0000000000 (ambient)**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Cellulose-2, 4 x 3.0mm, 10/Pk Part No.: AJ0-8366

- Holder Part No.: KJ0-4282

17312

**Products used in this application:****ANALYTES:****1** Warfarin

Retention Time: 5.7 min

2 Warfarin

Retention Time: 7.3 min

