

## Mianserin on Lux Cellulose-2

**Column:** Lux® 5 µm Cellulose-2, LC Column 250 x 4.6 mm, Ea

**Dimensions:** 250 x 4.6 mm ID

**Order No:** 00G-4457-E0

**Elution Type:** Isocratic

**Eluent A:** 0.1 % DEA in CH<sub>3</sub>CN

**Eluent B:** 0.1% DEA in IPA

| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B |
|-------------------|----------|------------|-------|-------|
| 1                 | 1        | 10         | 95    | 5     |

**Flow Rate:** 1 mL/min

**Col. Temp.:** ambient

**Detection:** Refractive Index @ 0.0000000000 (ambient)

**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.


- SecurityGuard Cartridges, Lux Cellulose-2, 4 x 3.0mm, 10/Pk Part No.: AJ0-8366

- Holder Part No.: KJ0-4282

17433



Products used in this application:

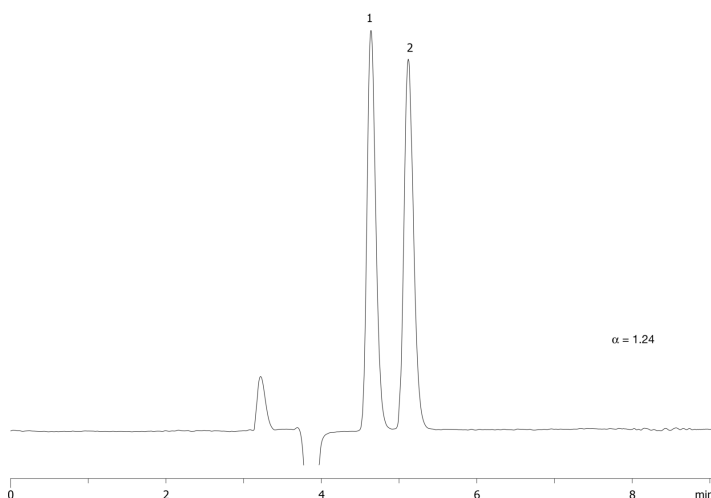

### ANALYTES:

**1** Mianserin

Retention Time: 4.6 min

**2** Mianserin

Retention Time: 5.1 min

