


**EPA 8270 Internal Standard Mix on ZB-MR2****Column:** Zebtron<sup>™</sup> ZB-MultiResidue<sup>™</sup>-2, GC Cap. Column 30 m x 0.25 mm x 0.20  $\mu$ m, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.25 mm x 0.2  $\mu$ m**Order No:** 7HG-G017-10**Oven Profile:****Carrier Gas:** Constant Flow Helium, 0.8 mL/min**Injection:** On-Column 0.000000000:1 0  $\mu$ L @ 43°C**Detection:** Refractive Index (335°C)**Analyst Note:** Analytes 10 ppm in Acetone

17596

**Products used in this application:****ANALYTES:**

- 1 1,4-Dichlorobenzene-d4
- 2 Naphthalene-d8
- 3 Acenaphthene-d10
- 4 Phenanthrene-d10
- 5 Chrysene-d12
- 6 Perylene-d12

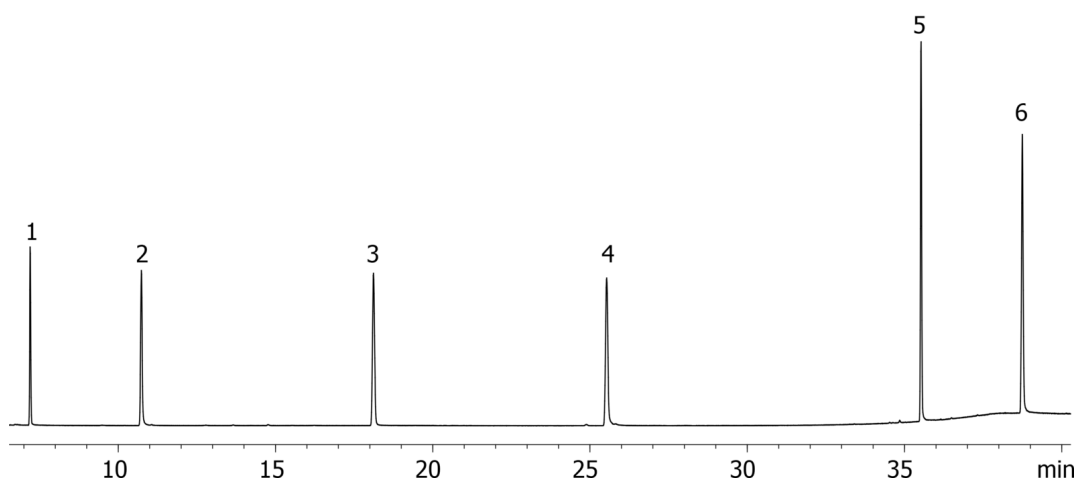