


Amphetamines by MSD on a ZB-1MS**Column:** Zebtron™ ZB-1MS, GC Cap. Column 10 m x 0.18 mm x 0.18 µm, Ea**Phase:** 100% Dimethylpolysiloxane**Dimensions:** 10 meters x 0.18 mm x 0.18 µm**Order No:** 7CD-G011-08**Oven Profile:****Carrier Gas:** Constant Flow Helium, 55 cm/sec**Injection:** Split 10.000000000:1 0 µL @ 250°C**Detection:** Refractive Index (320°C)**Analyst Note:** Analytes are 40ppm derivatized with HFBA
Sim ions:
(Quant ion*) D11-Amphetamine (*244,128): Amphetamine (*240,118,117): D14-Methamphetamine (*261,213):
Methamphetamine (*254,210,118)**Products used in this application:**

17617

**ANALYTES:**

- 1 Amphetamine-d11
- 2 Amphetamine
- 3 Methamphetamine-d14
- 4 Methamphetamine

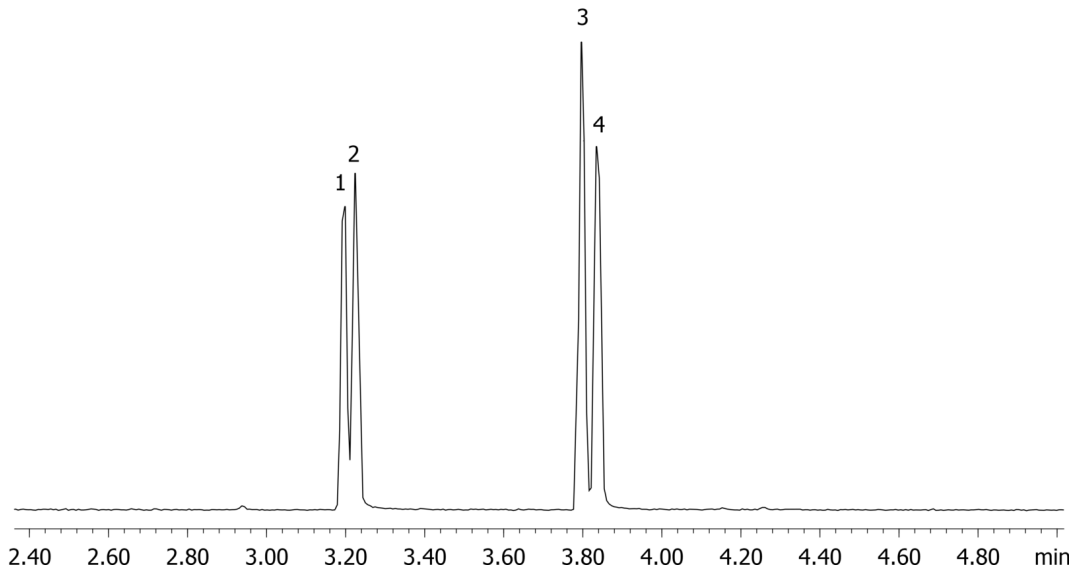