

Expanded amphetamines by MSD on a ZB-1MS**Column:** Zebron[™] ZB-1MS, GC Cap. Column 10 m x 0.18 mm x 0.18 µm, Ea**Phase:** 100% Dimethylpolysiloxane**Dimensions:** 10 meters x 0.18 mm x 0.18 µm**Order No:** 7CD-G011-08**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1 mL/min**Injection:** Split 50.0000000000:1 0 µL @ 280°C**Detection:** Refractive Index (230°C)**Analyst Note:** Analytes are 40ppm derivatized with HFBA


Sim ions:

(Quant ion*) D11-Amphetamine (*244,128): Amphetamine (*240,118,117): D14-Methamphetamine (*261,213):

Methamphetamine (*254,210,118): D5-MDA (*167,380,244): MDA (*162,375,240): D5-MDMA (*258,213,164):

**Products used in this application:**

17618

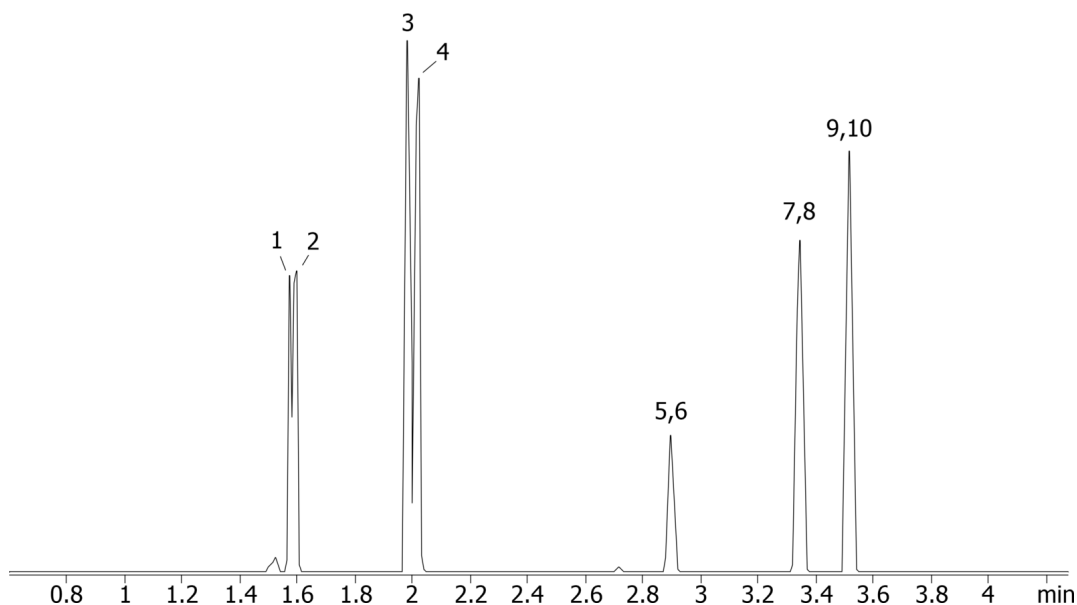

Expanded amphetamines by MSD on a ZB-1MS

ANALYTES:

- 1 Amphetamine-d11
- 2 Amphetamine
- 3 Methamphetamine-d14
- 4 Methamphetamine
- 5 MDA-d5
- 6 MDA
- 7 MDMA-d5
- 8 MDMA
- 9 MDEA-d5
- 10 MDEA

