


Benzene and butanol by MSD on a ZB-5MSi (d)**Column:** Zebron[™] ZB-5MSi, GC Cap. Column 30 m x 0.25 mm x 0.50 μ m, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.5 μ m**Order No:** 7HG-G018-17**Oven Profile:****Carrier Gas:** Constant Flow Helium, 55 cm/sec**Injection:** Split 50.0000000000:1 0 μ L @ 250°C**Detection:** Refractive Index (250°C)**Analyst Note:** Analytes are 50ppm in DMSO

17687

**Products used in this application:****ANALYTES:**

- 1 Benzene
- 2 Butanol

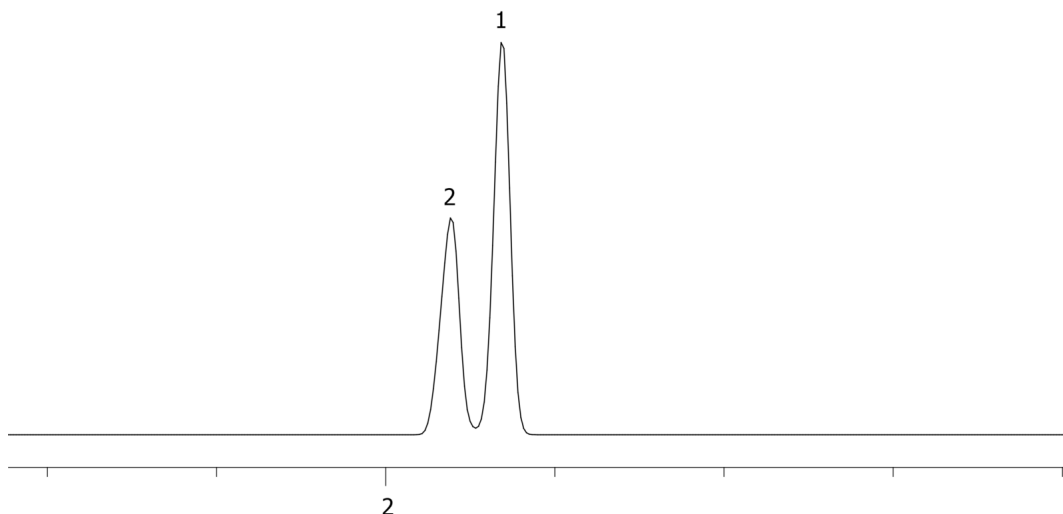