


Expanded amphetamines by MSD on a ZB-5MSi**Column:** Zebtron™ ZB-5MSi, GC Cap. Column 10 m x 0.18 mm x 0.18 µm, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 10 meters x 0.18 mm x 0.18 µm**Order No:** 7CD-G018-08**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1 mL/min**Injection:** Split 50.0000000000:1 0 µL @ 280°C**Detection:** Refractive Index (320°C)**Analyst Note:** Analytes are 40ppm derivatized with HFBA

Sim ions:

(Quant ion*) D11-Amphetamine (*244,128): Amphetamine (*240,118,117): D14-Methamphetamine (*261,213):
Methamphetamine (*254,210,118): D5-MDA (*167,380,244): MDA (*162,375,240): D5-MDMA (*258,213,164):**Products used in this application:**

17786

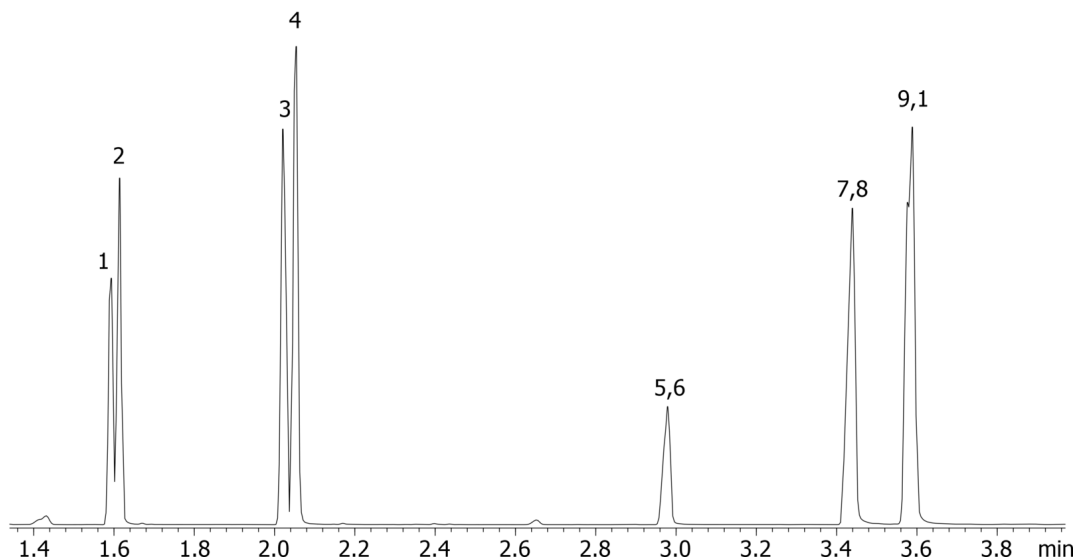

Expanded amphetamines by MSD on a ZB-5MSi

ANALYTES:

- 1 Amphetamine-d11
- 2 Amphetamine
- 3 Methamphetamine-d14
- 4 Methamphetamine
- 5 MDA-d5
- 6 MDA
- 7 MDMA-d5
- 8 MDMA
- 9 MDEA-d5
- 10 MDEA

