


Organophosphorous Pesticides Mix 2 on ZB-50**Column:** Zebtron™ ZB-50, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** 50% Phenyl 50% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G004-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 0.7 mL/min**Injection:** On-Column 0.000000000:1 0 µL @ 43°C**Detection:** Refractive Index (335°C)**Analyst Note:** Analytes are 10 ppm in Acetone**Products used in this application:**

Organophosphorous Pesticides Mix 2 on ZB-50

ANALYTES:

- 1 Acephate
- 2 Omethoate
- 3 Cadusafos
- 4 Monocrotophos
- 5 Dimethoate
- 6 Tris(chloropropyl) phosphate
- 7 Tris(chloropropyl) phosphate
- 8 Tolclofos-methyl
- 9 Methyl parathion
- 10 Pirimiphos-methyl
- 11 Fenthion
- 12 Parathion
- 13 Quinalphos
- 14 Methidathion
- 15 Prothiophos
- 16 Tribufos
- 17 Triazophos
- 18 EPN
- 19 Phosalone

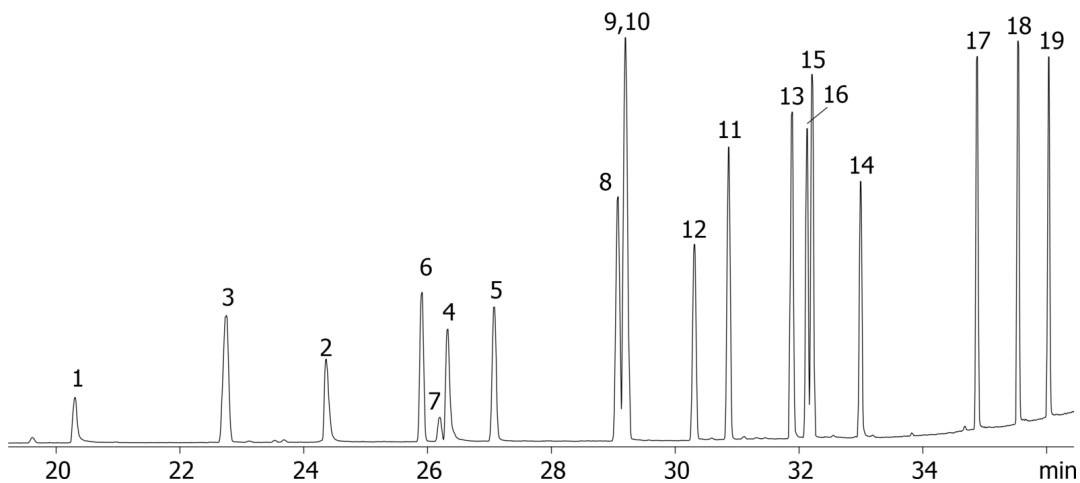