


EPA 8270 Internal Standard Mix on ZB-50**Column:** ZebtronTM ZB-50, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** 50% Phenyl 50% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G004-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 0.7 mL/min**Injection:** On-Column 0.000000000:1 0 μ L @ 43°C**Detection:** Refractive Index (335°C)**Analyst Note:** Analytes are 10 ppm in Acetone**Products used in this application:****ANALYTES:**

- 1 1,4-Dichlorobenzene-d4
- 2 Naphthalene-d8
- 3 Acenaphthene-d10
- 4 Phenanthrene-d10
- 5 Chrysene-d12
- 6 Perylene-d12

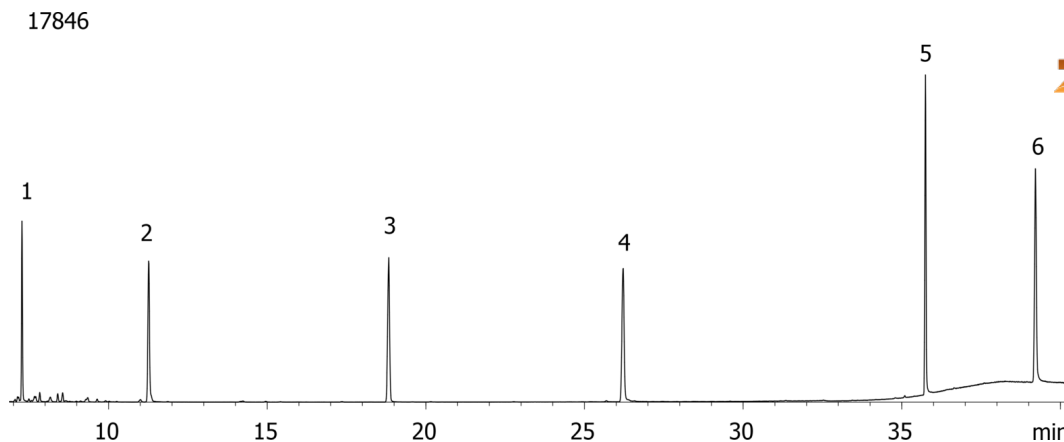