


Blood Alcohols on ZB-BAC-2**Column:** Zebtron[™] ZB-BAC-2, GC Cap. Column 30 m x 0.32 mm x 1.20 μ m, Ea**Phase:** Proprietary Blood Alcohol Phase**Dimensions:** 30 meters x 0.32 mm x 1.2 μ m**Order No:** 7HM-G022-25**Oven Profile:****Carrier Gas:** Constant Flow Helium, 80 cm/sec**Injection:** Split 5.000000000:1 0 μ L @ 140°C**Detection:** Refractive Index (220°C)**Analyst Note:** 0.050% Blood Alcohols Components in Water (1mL in 20mL Headspace Vial)

17898

**Products used in this application:****ANALYTES:**

- 1 Methanol
- 2 Acetaldehyde
- 3 Ethanol
- 4 Isopropanol
- 5 Acetone
- 6 t-Butanol
- 7 2-Butanol

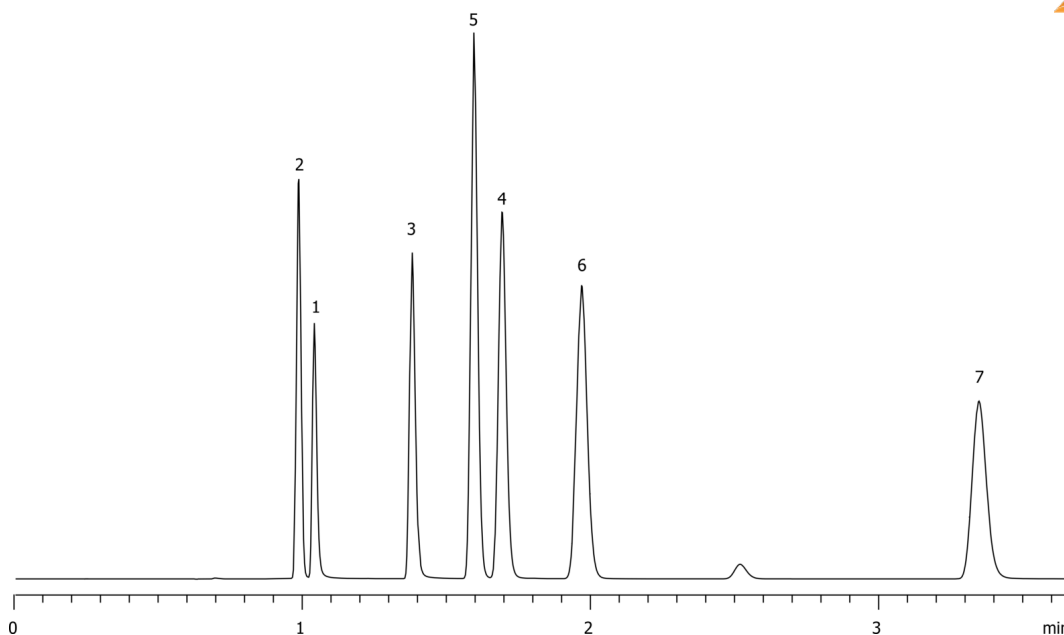