

ABN Mix on Synergi 4u Polar-RP 150x4.6


Column: SynergiTM 4 μ m Polar-RP 80 Å, LC Column 150 x 4.6 mm, Ea**Dimensions:** 150 x 4.6 mm ID**Order No:** 00F-4336-E0**Elution Type:** Gradient**Eluent A:** Phosphate Buffer, pH 2.5**Eluent B:** Acetonitrile**Eluent C:** 0**Eluent D:** 0

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B	Pct C	Pct D
	1	0	95	5	0	0
	2	8	20	80	0	0
	3	10	20	80	0	0
	4	10.5	95	5	0	0
	5	12	95	5	0	0

Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** Refractive Index @ 0.0000000000 (ambient)**Analyst Note:** SecurityGuardTM Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Polar-RP 4 x 3.0mm, 10/Pk Part No.: AJ0-6076

- Holder Part No.: KJ0-4282

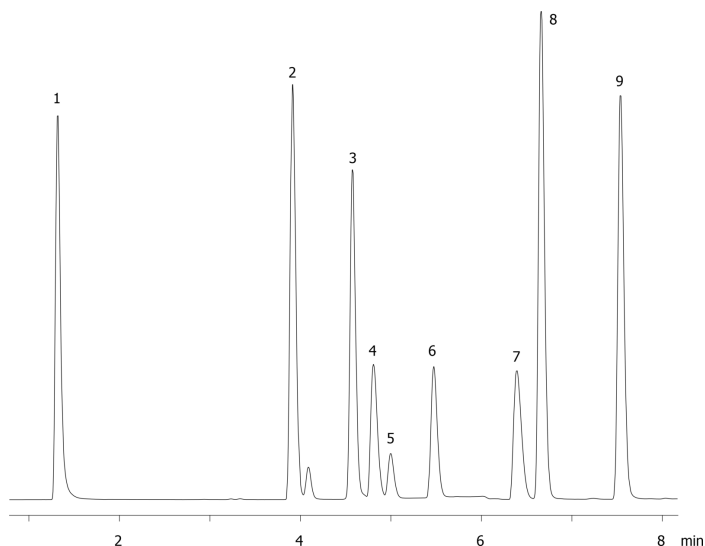

ANALYTES:

- 1 Pyridine
Retention Time: 1.3 min
- 2 Quinidine
Retention Time: 3.9 min
- 3 Sulfathiazole
Retention Time: 4.6 min
- 4 Triprolidine
Retention Time: 4.8 min
- 5 Benzyl alcohol
Retention Time: 5 min
- 6 Nortryptiline
Retention Time: 5.5 min
- 7 Phenol
Retention Time: 5.7 min

©2025 Phenomenex, Inc. All rights reserved.

For more information contact your Phenomenex Representative at support@phenomenex.com

 3-Methyl 4-benzoic acid
Phenomenex products are available worldwide.
Retention Time: 6.7 min
www.phenomenex.com

- 9 3-Methyl salicylaldehyde
Retention Time: 7.5 min

support@phenomenex.com