


**Blood Alcohols on ZB-BAC-1****Column:** Zebron<sup>™</sup> ZB-BAC-1, GC Cap. Column 30 m x 0.53 mm x 3.00  $\mu$ m, Ea**Phase:** Proprietary Blood Alcohol Phase**Dimensions:** 30 meters x 0.53 mm x 3  $\mu$ m**Order No:** 7HK-G021-36**Oven Profile:****Carrier Gas:** Constant Flow Helium, 80 cm/sec**Injection:** Split 0.8000000000:1 0  $\mu$ L @ 150°C**Detection:** Refractive Index (250°C)**Analyst Note:** Analytes 0.025% and internal standards 0.100% in water**Products used in this application:**

17982

**ANALYTES:**

- 1 Methanol
- 2 Acetaldehyde
- 3 Ethanol
- 4 Isopropanol
- 5 Acetone
- 6 t-Butanol
- 7 n-Propanol
- 8 2-Butanol

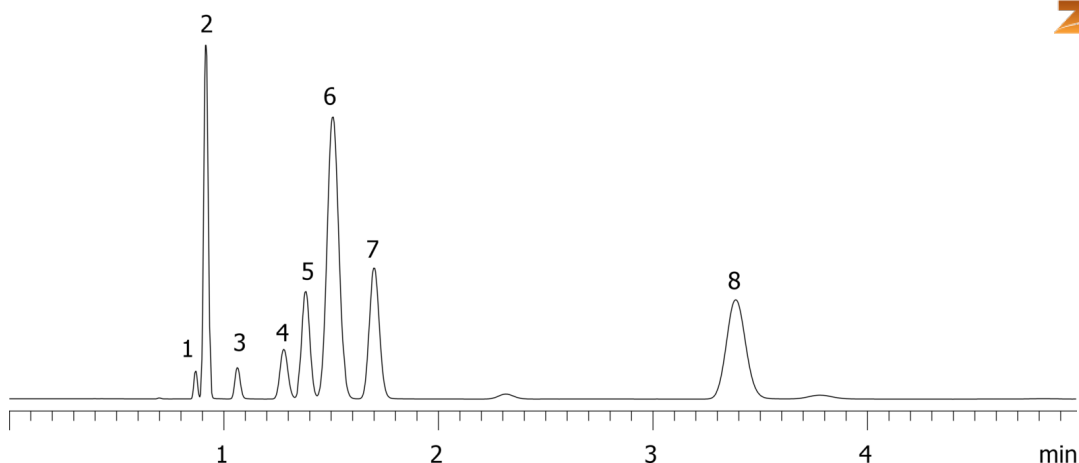