


Methadone from urine by Screen C on a ZB-Drug1**Column:** Zebtron™ ZB-Drug-1, GC Cap. Column 10 m x 0.18 mm x 0.18 µm, Ea**Phase:** Proprietary Drug Phase**Dimensions:** 10 meters x 0.18 mm x 0.18 µm**Order No:** 7CD-G023-08**Oven Profile:****Carrier Gas:** Constant Flow Helium, 55 cm/sec**Injection:** Split 20.0000000000:1 0 µL @ 280°C**Detection:** Refractive Index (320°C)**Analyst Note:** Analytes are 300ng/mL in urine
Sim ions:
(Quant ion*) D9-Methadone (*78,303): Methadone (*72,294,223)**Products used in this application:****ANALYTES:**

- 1 Methadone-d9
- 2 Methadone

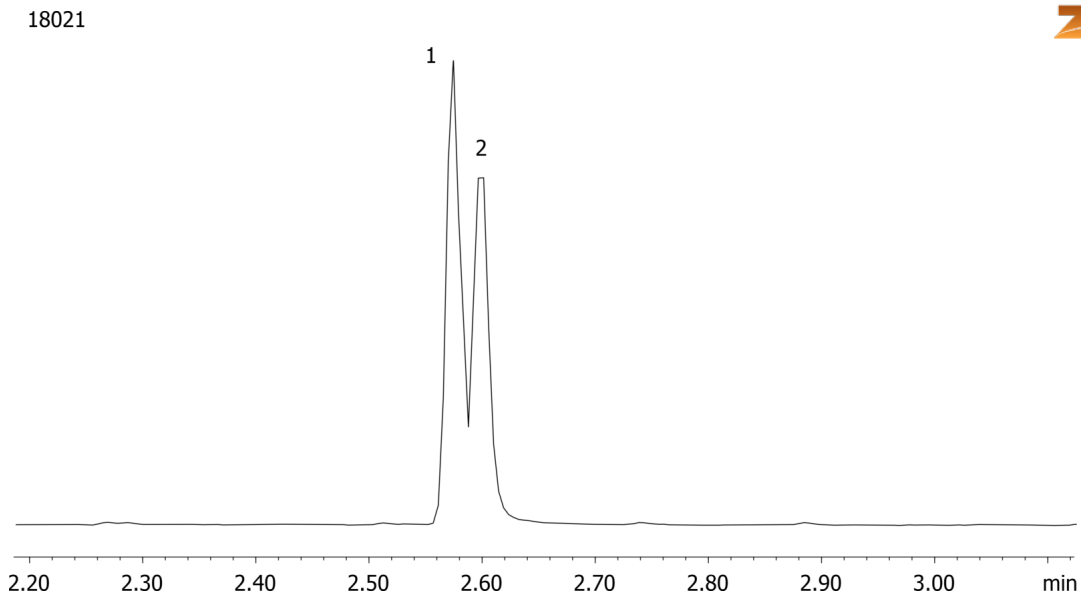