


6-MAM from urine by Screen C on a ZB-Drug1**Column:** Zebtron[™] ZB-Drug-1, GC Cap. Column 10 m x 0.18 mm x 0.18 µm, Ea**Phase:** Proprietary Drug Phase**Dimensions:** 10 meters x 0.18 mm x 0.18 µm**Order No:** 7CD-G023-08**Oven Profile:****Carrier Gas:** Constant Flow Helium, 0.7 mL/min**Injection:** Split 10.0000000000:1 0 µL @ 260°C**Detection:** Refractive Index (300°C)**Analyst Note:** Analytes are 10ng/mL in urine derivatized with BSTFA
Sim ions:
(Quant ion*) D3-6-MAM (*402,403): 6-MAM (*399,400,340)**Products used in this application:**

18022

**ANALYTES:**

- 1 6-MAM-d3
- 2 6-MAM

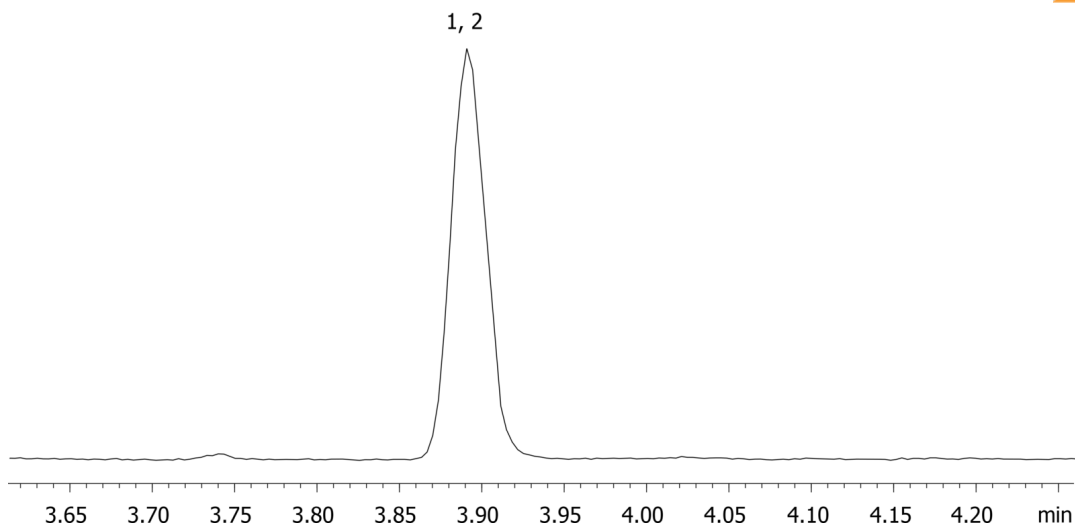