


Blood Alcohols using Hydrogen on ZB-BAC-2**Column:** Zebtron™ ZB-BAC-2, GC Cap. Column 30 m x 0.53 mm x 2.00 µm, Ea**Phase:** Proprietary Blood Alcohol Phase**Dimensions:** 30 meters x 0.53 mm x 2 µm**Order No:** 7HK-G022-32**Oven Profile:****Carrier Gas:** Constant Flow Hydrogen, 80 cm/sec**Injection:** Split 5.000000000:1 0 µL @ 150°C**Detection:** Refractive Index (250°C)**Analyst Note:** All analytes are 0.100% in water**Products used in this application:****ANALYTES:**

- 1 Methanol
- 2 Acetaldehyde
- 3 Ethanol
- 4 Isopropanol
- 5 Acetone
- 6 n-Propanol

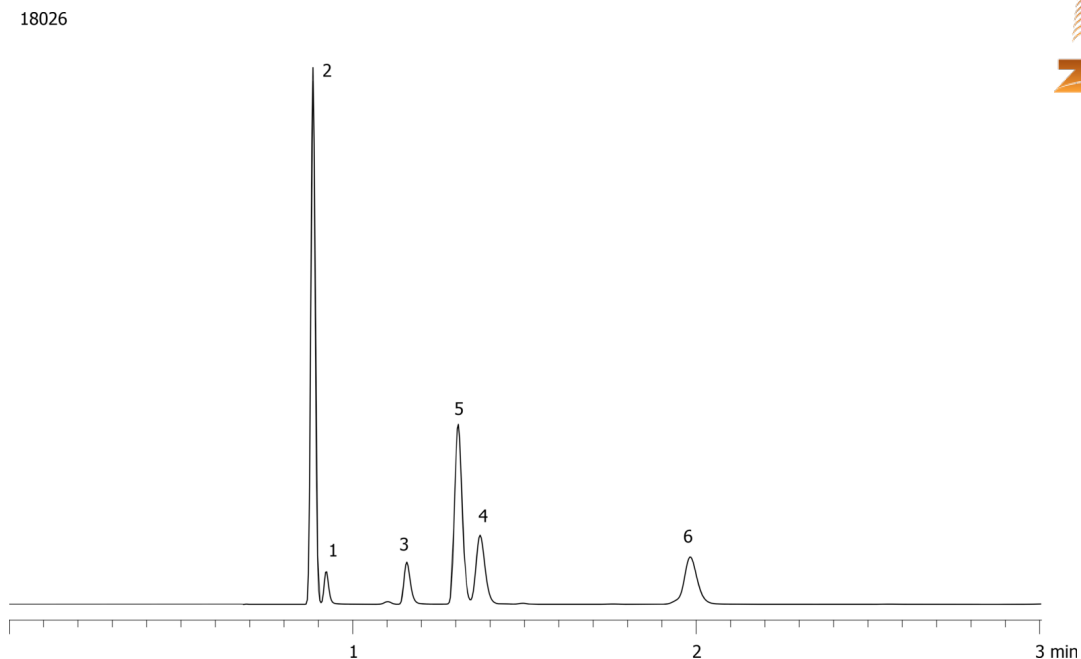