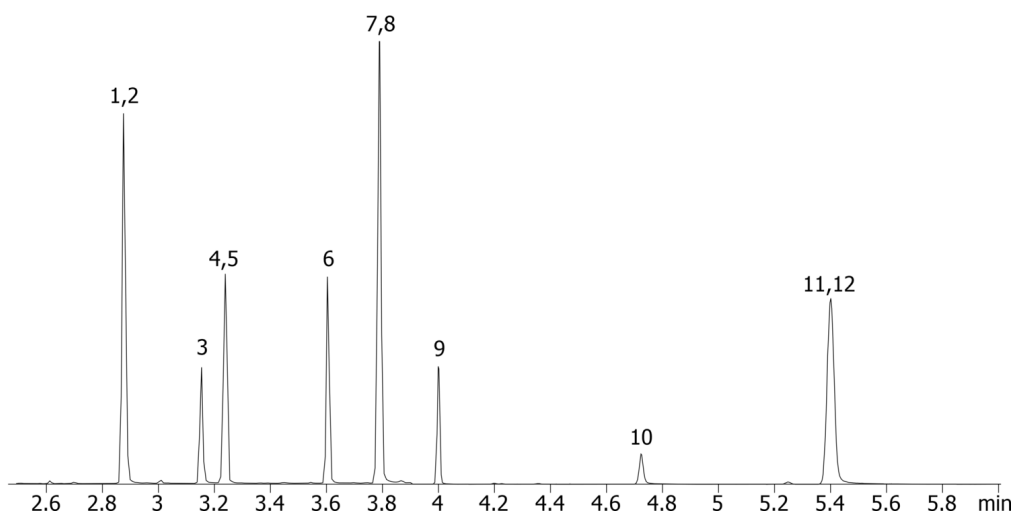


Benzodiazepines from urine by Screen C on a ZB-Drug1**Column:** Zebtron™ ZB-Drug-1, GC Cap. Column 10 m x 0.18 mm x 0.18 µm, Ea**Phase:** Proprietary Drug Phase**Dimensions:** 10 meters x 0.18 mm x 0.18 µm**Order No:** 7CD-G023-08**Oven Profile:****Carrier Gas:** Constant Flow Helium, 0.7 mL/min**Injection:** Split 10.0000000000:1 0 µL @ 260°C**Detection:** Refractive Index (320°C)**Analyst Note:** Analytes are 300ng/mL in urine derivatized with MTBSTFA

Sim ions:

(Quant ion*) D5-Nordiazepam (*332,334): Nordiazepam (*327,328,383): D5-Oxazepam (*462,519): Oxazepam (*457,513,459): D5-Temazepam (*362,288,390): Temazepam (*357,359,385): Lorazepam (*491,513,493):

18055

**Products used in this application:**

Benzodiazepines from urine by Screen C on a ZB-Drug1

ANALYTES:

- 1 Nordiazepam-d5
- 2 Nordiazepam
- 3 Diazepam
- 4 Oxazepam-d5
- 5 Oxazepam
- 6 Lorazepam
- 7 Temazepam-d5
- 8 Temazepam
- 9 Clonazepam
- 10 Alprazolam
- 11 alpha-Hydroxyalprazolam-d5
- 12 alpha-Hydroxyalprazolam

