

HPLC Application

ID No.: 18099

Narcotics 8-mix on Luna 5 μ m C18(2) 150 x 4.6

Column: Luna[®] 5 μ m C18(2) 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4252-E0

Elution Type: Gradient

Eluent A: 10mM NH₄OAc pH5.5

Eluent B: ACN

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	3	95	5
	3	23	60	40

Flow Rate: 1 mL/min

Col. Temp.: 45 °C

Detection: Refractive Index @ 0.0000000000 (ambient)

Analyst Note: SecurityGuard[™] Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm ID, 10/Pk Part No.: AJ0-4287


- Holder Part No.: KJ0-4282



Products used in this application:



18099

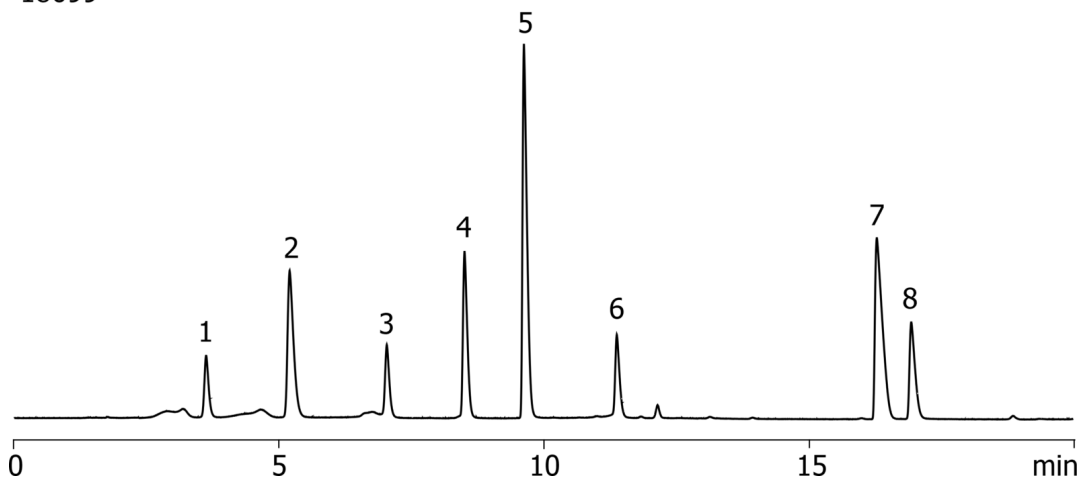

ANALYTES:

- 1 Normorphine
Retention Time: 3.648 min
- 2 Morphine
Retention Time: 5.222 min
- 3 Hydromorphone
Retention Time: 7.062 min
- 4 Norcodeine
Retention Time: 8.517 min
- 5 Codeine
Retention Time: 9.637 min
- 6 Hydrocodone
Retention Time: 11.386 min
- 7 Cocaine
Retention Time: 16.28 min
- 8 Norcocaine
Retention Time: 16.93 min

