

Methyl p-tolyl sulfoxide on Lux 5 μ m Amylose-2 250 x 4.6mm in NP**Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** 0.1 % DEA in Hexane**Eluent B:** 0.1% DEA in IPA

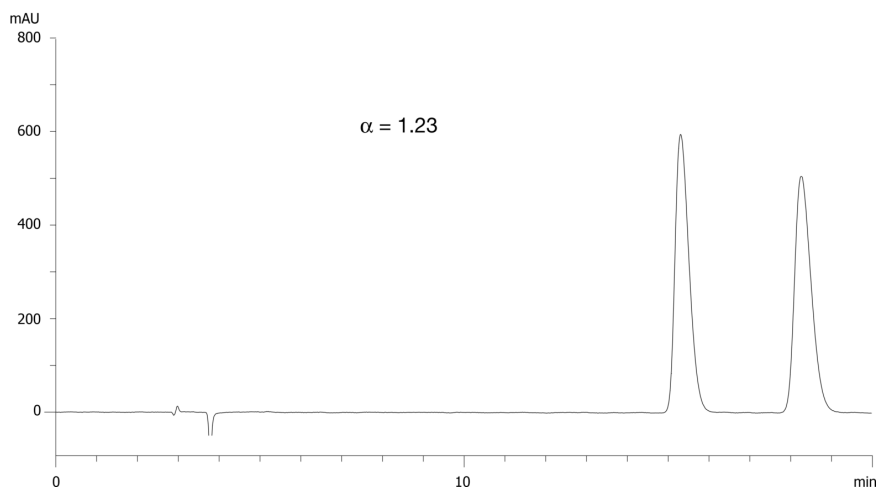
Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	10	80	20

Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** Refractive Index @ 0.0000000000 (ambient)**Analyst Note:** SecurityGuardTM Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Amylose-2, 4 x 3.0mm, 10/Pk Part No.: AJ0-8470

- Holder Part No.: KJ0-4282

18109

**ANALYTES:****1** Methyl p-tolyl sulfoxide

Retention Time: 15.3 min

2 Methyl p-tolyl sulfoxide

Retention Time: 18.26 min

