


Biodiesel Esters - PME on ZB-FFAP**Column:** Zebron[™] ZB-FFAP, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** Nitroterephthalic Acid Modified Polyethylene Glycol**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G009-11**Oven Profile:****Carrier Gas:** Constant Pressure Hydrogen, 1.6 mL/min**Injection:** Split 50.0000000000:1 0 μ L @ 250°C**Detection:** Refractive Index (270°C)**Analyst Note:** Courtesy of ASG mbH, Germany**Products used in this application:**

Biodiesel Esters - PME on ZB-FFAP

ANALYTES:

- 1 C8
- 2 C10
- 3 C12
- 4 C14
- 5 C16:0
- 6 C16:1
- 7 internal standard
- 8 C18:0
- 9 C18:1
- 10 C18:2
- 11 C18:3
- 12 C20:0
- 13 C20:1
- 14 C22:0
- 15 C23
- 16 C24:0

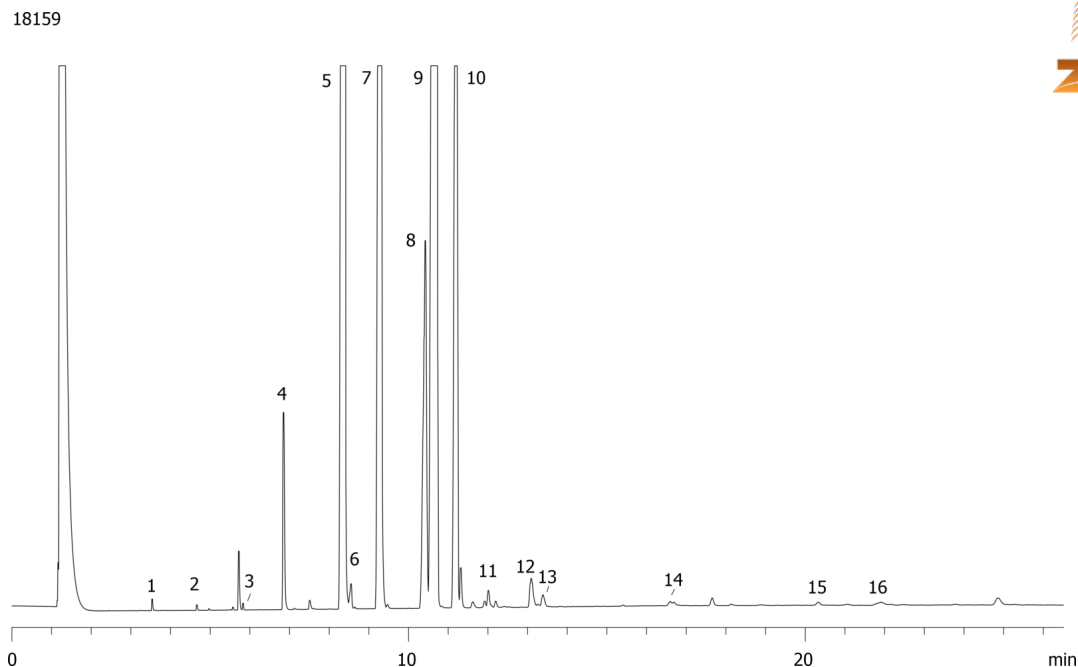