


Biodiesel Esters - RME on ZB-FFAP**Column:** Zebtron™ ZB-FFAP, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** Nitroterephthalic Acid Modified Polyethylene Glycol**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G009-11**Oven Profile:****Carrier Gas:** Constant Pressure Hydrogen, 1.6 mL/min**Injection:** Split 50.0000000000:1 0 µL @ 250°C**Detection:** Refractive Index (270°C)**Analyst Note:** Courtesy of ASG mbH, Germany**Products used in this application:**

Biodiesel Esters - RME on ZB-FFAP

ANALYTES:

- 1 C8
- 2 C12
- 3 C14
- 4 C16:0
- 5 C16:1
- 6 internal standard
- 7 C18:0
- 8 C18:1
- 9 C18:2
- 10 C18:3
- 11 C20:0
- 12 C20:1
- 13 C22:0
- 14 C23
- 15 C24:0

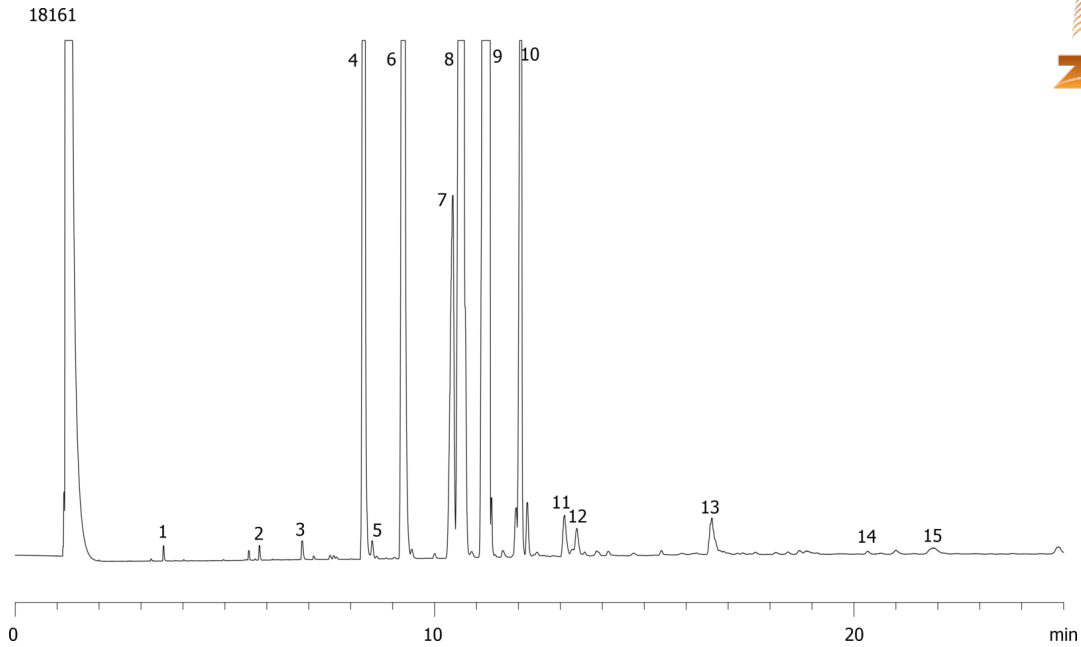