


European Multipesticide Analysis by MS on ZB-5**Column:** Zebron[™] ZB-5, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G002-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1 mL/min**Injection:** Splitless 0.000000000:1.0 µL @ 240°C**Detection:** Refractive Index (280°C)**Analyst Note:** Analytes are 1ppm. Pulsed Splitless injection @ 300kPa for 1.0 min.
Application courtesy of Frau Dr. F Klose at IVPT GmbH**Products used in this application:**

18182

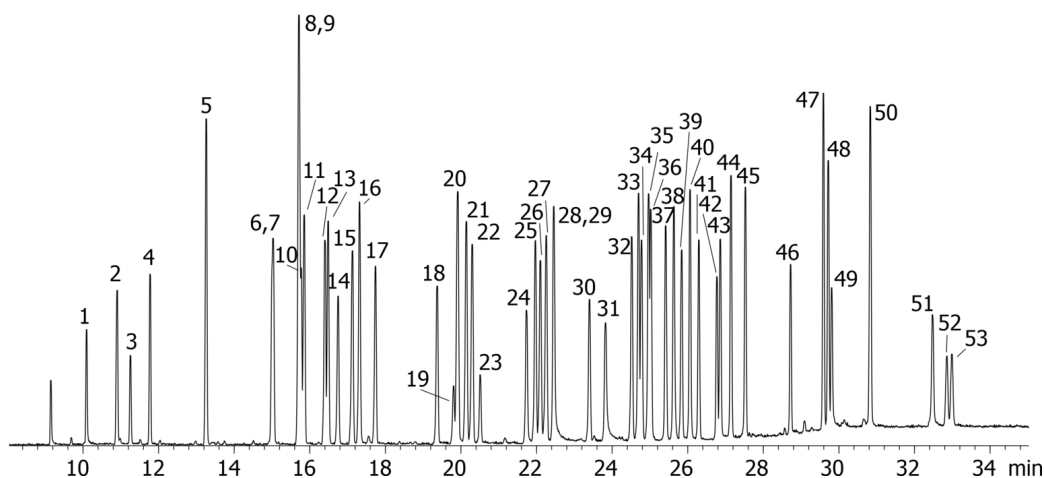

European Multipesticide Analysis by MS on ZB-5**ANALYTES:**

- | | |
|-----------------------------|----------------------------|
| 1 Dichlobenil | 41 Mefenpyr-diethyl |
| 2 Mevinphos | 42 EPN |
| 3 Nitrapyrin | 43 Picolinafen |
| 4 Methacrifos | 44 Fenamidone |
| 5 Propachlor | 45 Tetradifon |
| 6 Dicloran | 46 Pyrazophos |
| 7 Penntachloranisol | 47 Pyridaben |
| 8 Propetamphos | 48 Fluquinconazole |
| 9 Cyanophos | 49 Prochloraz |
| 10 Profluralin | 50 Halfenprox |
| 11 Fonofos | 51 Flumioxazin |
| 12 Tefluthrin | 52 tau-Fluvalinate |
| 13 Isazofos | 53 tau-Fluvalinate |
| 14 Iprobenfos | |
| 15 Pentachloranilin | |
| 16 Dichlofenthion | |
| 17 Transfluthrin | |
| 18 Flufenacet | |
| 19 Fosthiazate | |
| 20 Fosthiazate | |
| 21 Isofenphos-methyl | |
| 22 Metazachlor | |
| 23 Chlozolate | |
| 24 Climbazol | |
| 25 Chlorfenson | |
| 26 Hexaconazole | |
| 27 Prothiophos | |
| 28 Fludioxonil | |
| 29 Uniconazole | |
| 30 Chlorfenapyr | |
| 31 Diniconazole | |
| 32 Sulprofos | |
| 33 Isoxadifen-ethyl | |
| 34 Famphur | |
| 35 Cyanofenphos | |
| 36 Quinoxifen | |
| 37 Fluopicolide | |
| 38 Nuarimol | |
| 39 Diflufenican | |
| 40 Fluotrimazol | |

