


EPA Methods 508 and 8081A by MSD on ZB-35HT**Column:** Zebron[™] ZB-35HT, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** 35% Phenyl 65% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G025-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 0.5 mL/min**Injection:** Split 10.0000000000:1 0 µL @ 200°C**Detection:** Refractive Index (350°C)**Analyst Note:** All analytes were 10 ppm in iso-octane**Products used in this application:**

18199

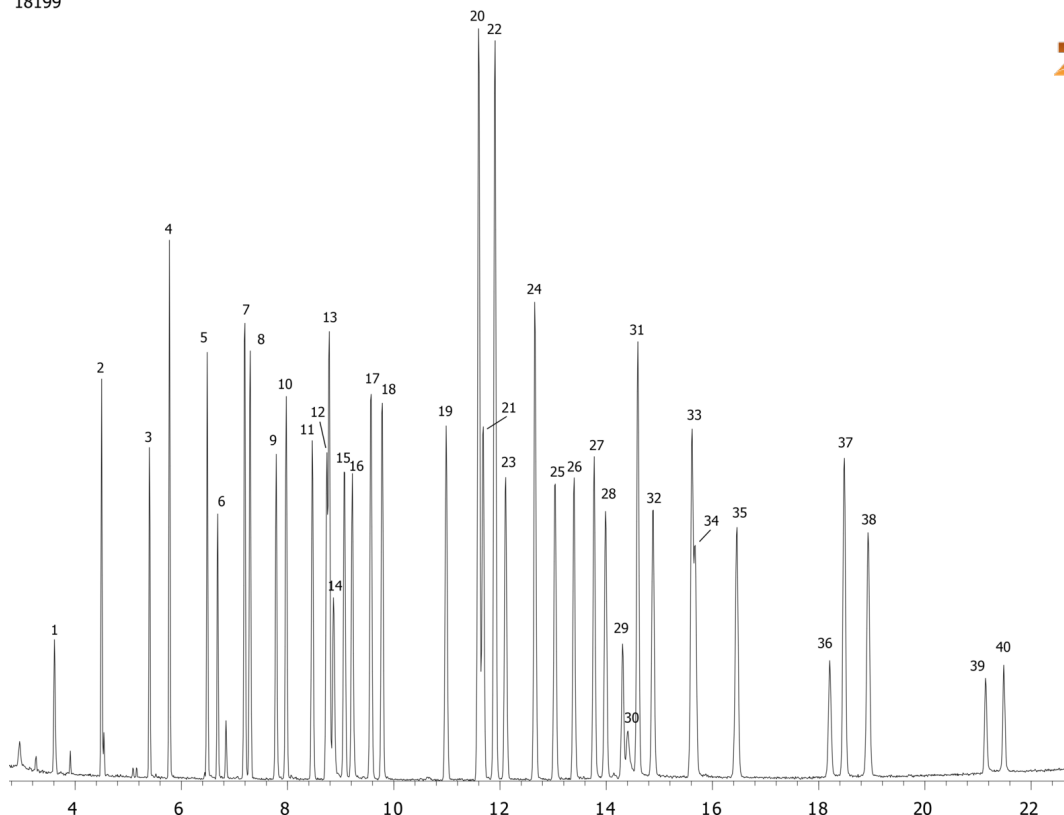

EPA Methods 508 and 8081A by MSD on ZB-35HT

ANALYTES:

- 1 DBCP
- 2 Hexachlorocyclopentadiene
- 3 Etridiazole
- 4 Chloroneb
- 5 Propachlor
- 6 Di-allate
- 7 Hexachlorobenzene
- 8 alpha-BHC
- 9 PCNB
- 10 gamma-BHC
- 11 b-BHC
- 12 Alachlor
- 13 Heptachlor
- 14 Dichlone
- 15 delta-BHC
- 16 Chlorothalonil
- 17 Aldrin
- 18 Dacthal
- 19 Heptachlor epoxide
- 20 trans-Chlordane
- 21 trans-Nonachlor
- 22 cis-Chlordane
- 23 Endosulfan I
- 24 4,4'-DDE
- 25 Dieldrin
- 26 Chloropropylate
- 27 Chlorobenzilate
- 28 Endrin
- 29 Nitrofen
- 30 Kepone
- 31 4,4'-DDD
- 32 Endosulfan II
- 33 4,4'-DDT
- 34 Endrin aldehyde
- 35 Endosulfan sulfate
- 36 Captafol
- 37 Methoxychlor
- 38 Endrin ketone
- 39 Permethrin
- 40 Permethrin

