

Simultaneous Analysis of Melamine and Cyanuric acid on Luna 3u 100x2.0 HILIC

Column: Luna® 3 µm HILIC 200 Å, LC Column 100 x 2 mm, Ea

Dimensions: 100 x 2 mm ID

Order No: 00D-4449-B0

Elution Type: Isocratic

Eluent A: 90/10- ACN/100 mM ammonium formate, pH 3.2

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100


Flow Rate: 0.4 mL/min

Col. Temp.: ambient

Detection: Refractive Index @ 0.0000000000 (ambient)

Analyst Note: LC/MS/MS Condition:

18289

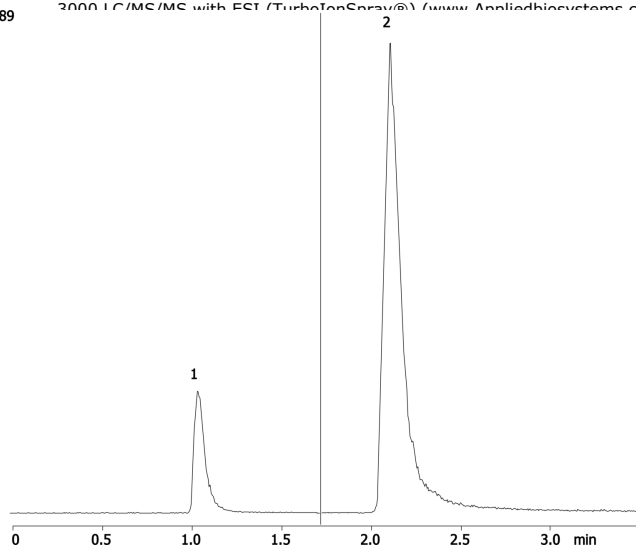

ANALYTES:

1 Cyanuric acid

Retention Time: 1.03 min

2 Melamine

Retention Time: 2.11 min



Products used in this application:

