


NSO (Nitrogen/Sulfur/Oxygen) Pesticides Mix on ZB-5MSi**Column:** Zebron[™] ZB-5MSi, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G018-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 0.8 mL/min**Injection:** On-Column 0.000000000:1 0 μ L @ 43°C**Detection:** Refractive Index (335°C)**Analyst Note:** Analytes 10 ppm in Acetone**Products used in this application:**

NSO (Nitrogen/Sulfur/Oxygen) Pesticides Mix on ZB-5MSi

ANALYTES:

- 1 Metaldehyde
- 2 Biphenyl
- 3 o-Phenylphenol
- 4 Diphenylamine
- 5 Ethoxyquin
- 6 Pyrimethanil
- 7 Pirimicarb
- 8 Metalaxyl
- 9 Carbaryl
- 10 Methiocarb
- 11 Diethofencarb
- 12 MGK-264
- 13 MGK-264
- 14 Cyprodinil
- 15 Thiabendazole
- 16 Isoprothiolane
- 17 Buprofezin
- 18 Bupirimate
- 19 Flusilazole
- 20 Diniconazole
- 21 Oxadixyl
- 22 Benalaxyl
- 23 Propargite
- 24 Piperonyl butoxide
- 25 Fenpropathrin
- 26 Pyriproxifen
- 27 Bitertanol
- 28 Bitertanol
- 29 Azoxystrobin

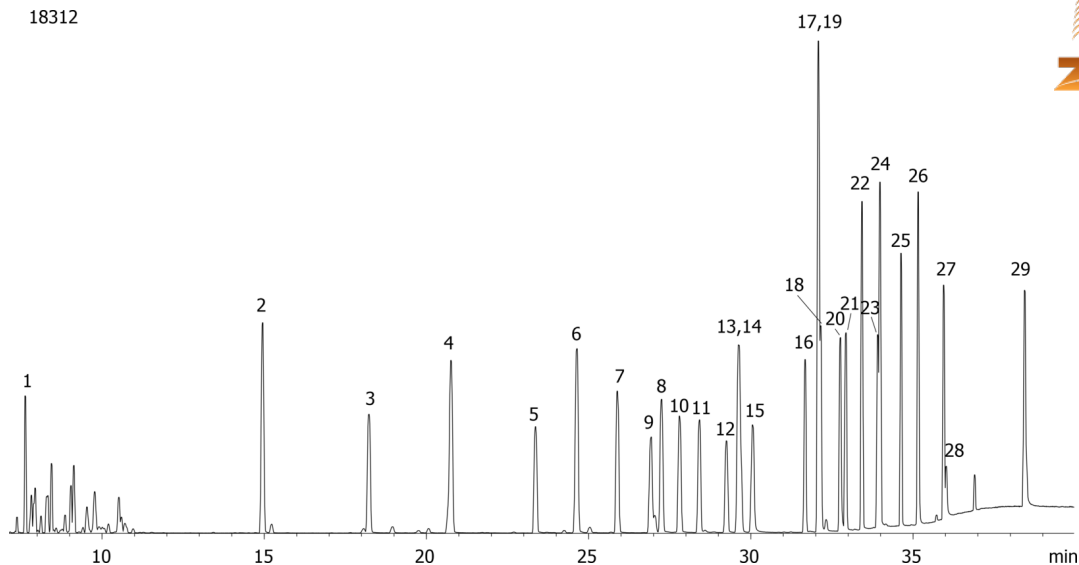