

# HPLC Application

ID No.: 18456

## Methadone on Lux 5 um Cellulose-2 250 x 4.6

**Column:** Lux® 5 µm Cellulose-2, LC Column 250 x 4.6 mm, Ea

**Dimensions:** 250 x 4.6 mm ID

**Order No:** 00G-4457-E0

**Elution Type:** Isocratic

**Eluent A:** 0.5 % DEA in Hexane

**Eluent B:** 0.5 % DEA in Ethanol

| Gradient | Step No. | Time (min) | Pct A | Pct B |
|----------|----------|------------|-------|-------|
| Profile: | 1        | 0          | 99    | 1     |

**Flow Rate:** 1 mL/min

**Col. Temp.:** 20 °C

**Detection:** Refractive Index @ 0.0000000000 (ambient)

**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Cellulose-2, 4 x 3.0mm, 10/Pk Part No.: AJ0-8366


- Holder Part No.: KJ0-4282



Products used in this application:



18456

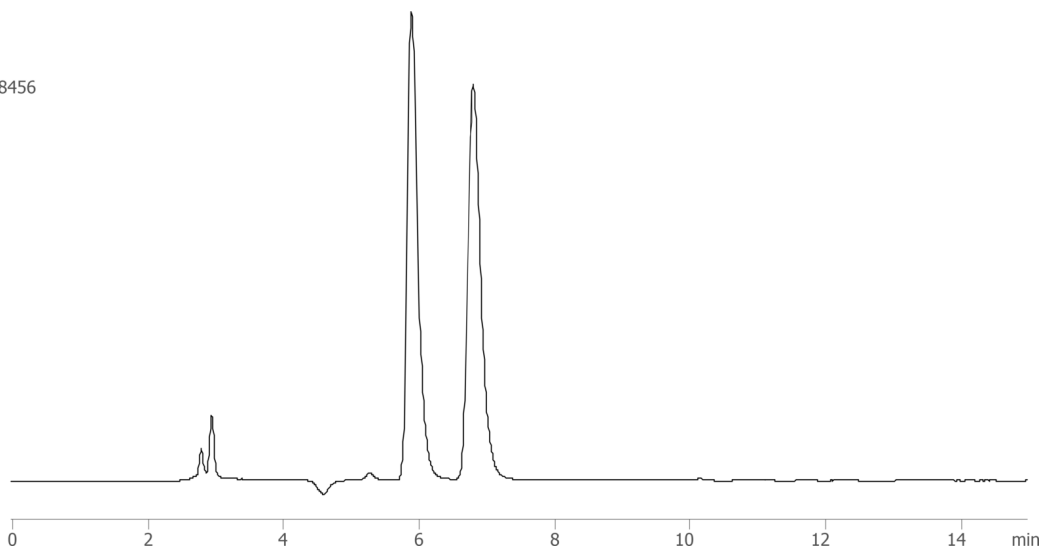

### ANALYTES:

1 Methadone

Retention Time: 5.892 min

2 Methadone

Retention Time: 6.802 min

