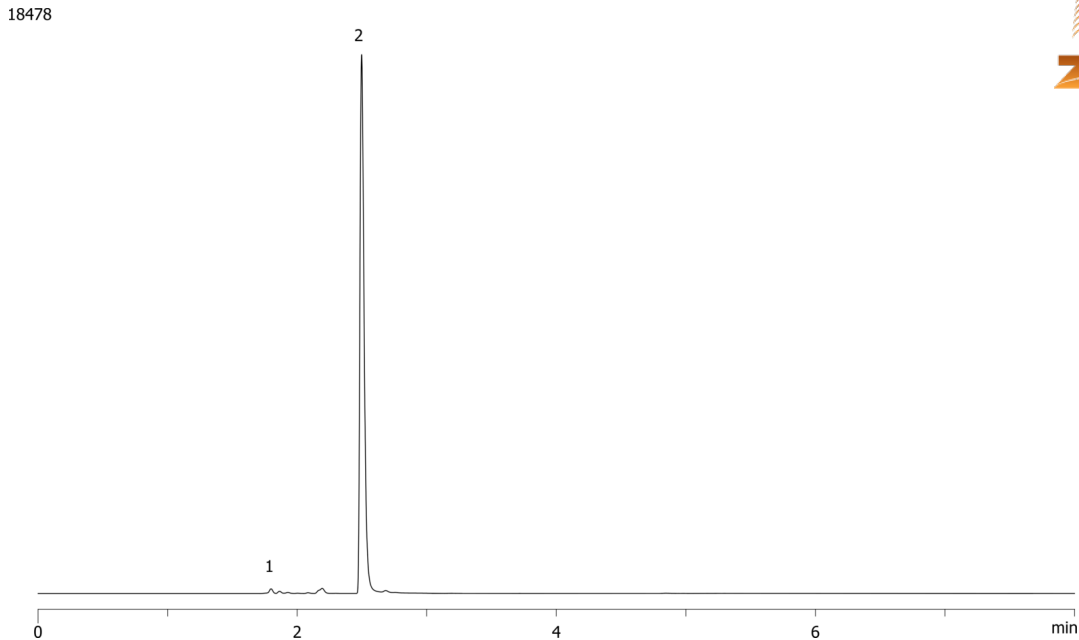


BioMethanol Analysis - RME on ZB-FFAP

Column: Zebron™ ZB-FFAP, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea
Phase: Nitroterephthalic Acid Modified Polyethylene Glycol
Dimensions: 30 meters x 0.25 mm x 0.25 µm
Order No: 7HG-G009-11
Oven Profile:

Carrier Gas: Constant Flow Hydrogen, 1 mL/min
Injection: Split 50.0000000000:1 0 µL @ 250°C
Detection: Refractive Index (300°C)
Analyst Note: Courtesy of ASG (www.asg-analytik.de), a specialized Biofuel Laboratory!

**Products used in this application:****ANALYTES:**

- 1 Hexaninoel
- 2 Methanol-FAME

