


BioMethanol Analysis - SME on ZB-FFAP**Column:** Zebron[™] ZB-FFAP, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** Nitroterephthalic Acid Modified Polyethylene Glycol**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G009-11**Oven Profile:****Carrier Gas:** Constant Flow Hydrogen, 1 mL/min**Injection:** Split 50.0000000000:1 0 µL @ 250°C**Detection:** Refractive Index (300°C)**Analyst Note:** Courtesy of ASG (www.asg-analytik.de), a specialized Biofuel Laboratory!**Products used in this application:**

18479

**ANALYTES:**

- 1 Hexaninoel
- 2 Methanol-FAME

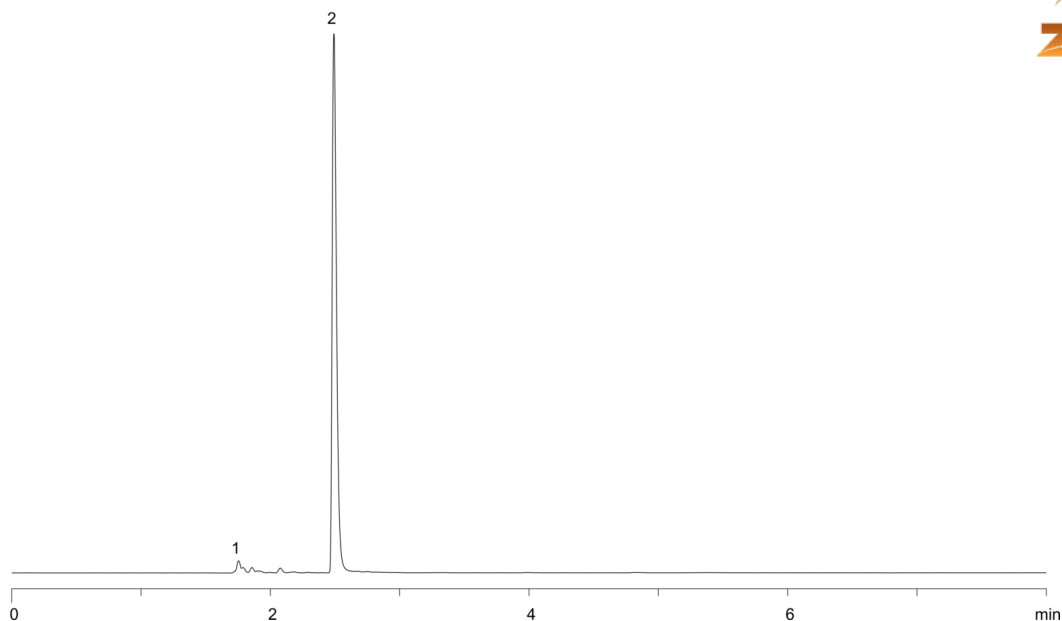