

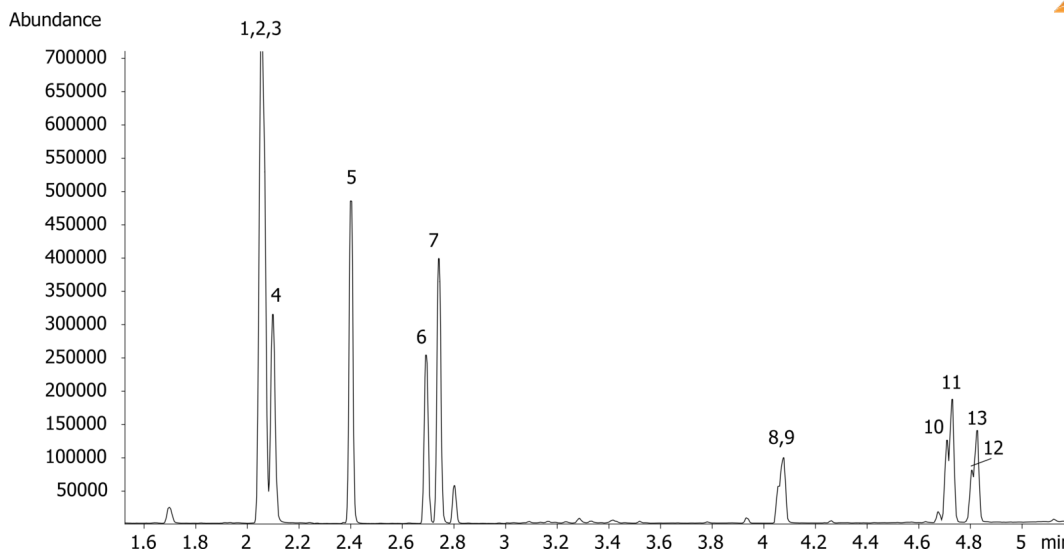
Alternative amphetamines by MSD on a ZB-Drug1**Column:** Zebron[™] ZB-Drug-1, GC Cap. Column 10 m x 0.18 mm x 0.18 µm, Ea**Phase:** Proprietary Drug Phase**Dimensions:** 10 meters x 0.18 mm x 0.18 µm**Order No:** 7CD-G023-08**Oven Profile:****Carrier Gas:** Constant Flow Helium, 0.7 mL/min**Injection:** Split 10.0000000000:1 0 µL @ 260°C**Detection:** Refractive Index (320°C)**Analyst Note:** Analytes are 20ppm derivatized with HFBA

Sim ions:

(Quant ion*) D11-Amphetamine (*244,128): Amphetamine (*240,118,117): D14-Methamphetamine (*261,213):

Methamphetamine (*254,210,118): D5-MDA (*167,380,244): MDA (*162,375,240): D5-MDMA (*258,213,164):

MDMA (*254,210,389): D5-MDEA (*273,241,408): MDEA (*268,240,403): Phentermine (*254,214,132):

**Products used in this application:**

Alternative amphetamines by MSD on a ZB-Drug1

ANALYTES:

- 1** Phenylpropanolamine
- 2** Phentermine
- 3** Amphetamine-d11
- 4** Amphetamine
- 5** Ephedrine
- 6** Methamphetamine-d14
- 7** Methamphetamine
- 8** MDA-d5
- 9** MDA
- 10** MDMA-d5
- 11** MDMA
- 12** MDEA-d5
- 13** MDEA

