


Polycyclic Aromatic Hydrocarbons (PAHs) in Water by Zebron ZB-35**Column:** Zebron[™] ZB-35, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** 35% Phenyl 65% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G003-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1.2 mL/min**Injection:** On-Column 0.000000000:1 0 µL @ 83°C**Detection:** Refractive Index (360°C)**Analyst Note:** Analytes are 10 ppm in dichloromethane**Products used in this application:**

Polycyclic Aromatic Hydrocarbons (PAHs) in Water by Zebron ZB-35**ANALYTES:**

- 1 Naphthalene
- 2 2-Methylnaphthalene
- 3 Acenaphthylene
- 4 Acenaphthene
- 5 Fluorene
- 6 Phenanthrene
- 7 Anthracene
- 8 Fluoranthene
- 9 Pyrene
- 10 Benzo[a]anthracene
- 11 Chrysene
- 12 Benzo[b]fluoranthene
- 13 Benzo[k]fluoranthene
- 14 Benzo[j]fluoranthene
- 15 Benzo[a]pyrene
- 16 3-Methylcholanthrene
- 17 Dibenz[a,h]acridine
- 18 Dibenz[a,,j]acridine
- 19 Indeno[1,2,3-cd]pyrene
- 20 Dibenz[a,h]anthracene
- 21 Benzo[g,h,i]perylene
- 22 7H-Dibenzo[c,g]carbazole
- 23 Dibenzo[a,e]pyrene
- 24 Dibenzo[a,i]pyrene
- 25 Dibenzo[a,h]pyrene

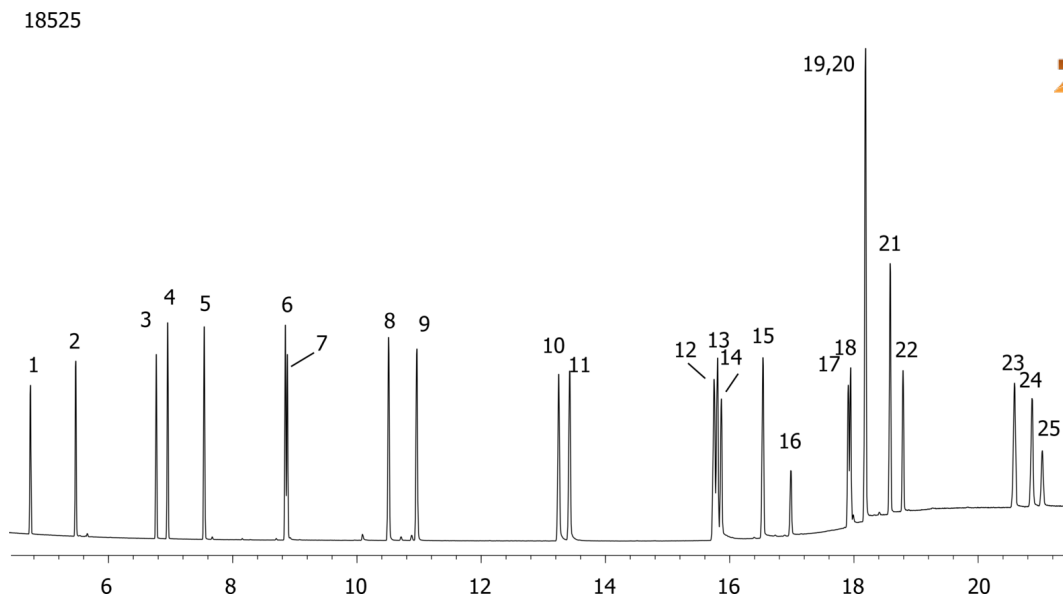