

FMOC-DL-Isoleucine-(boc)-OH on Lux 5 µm Cellulose-2 250 x 4.6mm in NP**Column:** Lux® 5 µm Cellulose-2, LC Column 250 x 4.6 mm, Ea**Dimensions:** 250 x 4.6 mm ID**Order No:** 00G-4457-E0**Elution Type:** Isocratic**Eluent A:** 0.1 % Formic acid in Hexane**Eluent B:** 0.1% Formic acid in IPA

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	15	80	20

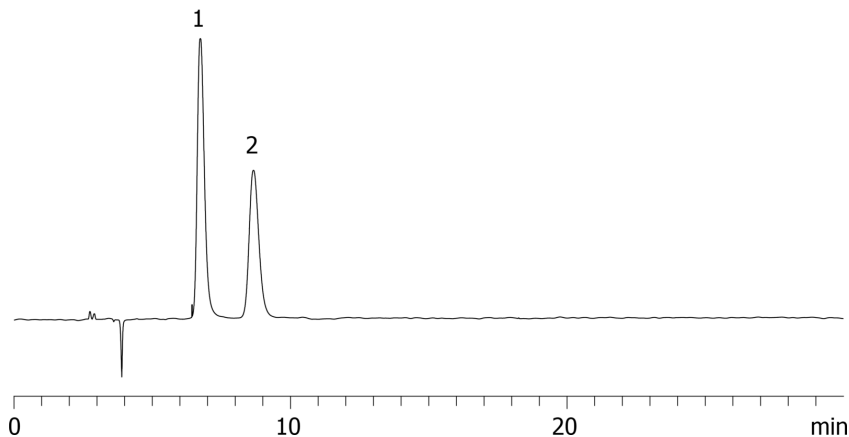
Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** Refractive Index @ 0.0000000000 (ambient)**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Cellulose-2, 4 x 3.0mm ID, 10/Pk Part No.: AJ0-8366

- Holder Part No.: KJ0-4282

**Products used in this application:**

18727

**ANALYTES:**

- 1** FMOC-Ile-OH
Retention Time: 6.7 min
- 2** FMOC-Ile-OH
Retention Time: 8.6 min

