

FMOC-DL-Serine-(boc)-OH on Lux 5 µm Cellulose-2 250 x 4.6mm in NP

Column: Lux® 5 µm Cellulose-2, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4457-E0

Elution Type: Isocratic

Eluent A: 0.1% TFA in Hexane

Eluent B: 0.1% TFA in IPA

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
1	1	20	85	15

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: Refractive Index @ 0.0000000000 (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.


- SecurityGuard Cartridges, Lux Cellulose-2, 4 x 3.0mm ID, 10/Pk Part No.: AJ0-8366

- Holder Part No.: KJ0-4282

18735



Products used in this application:

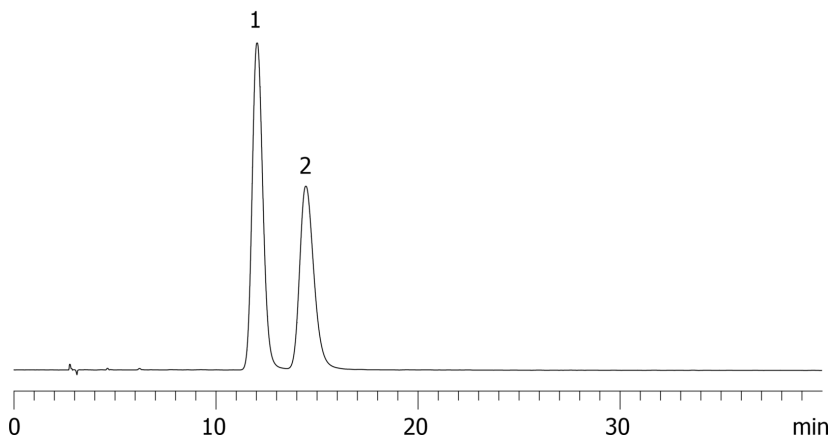

ANALYTES:

1 FMOC-Ser-OH

Retention Time: 12 min

2 FMOC-Ser-OH

Retention Time: 14.4 min

