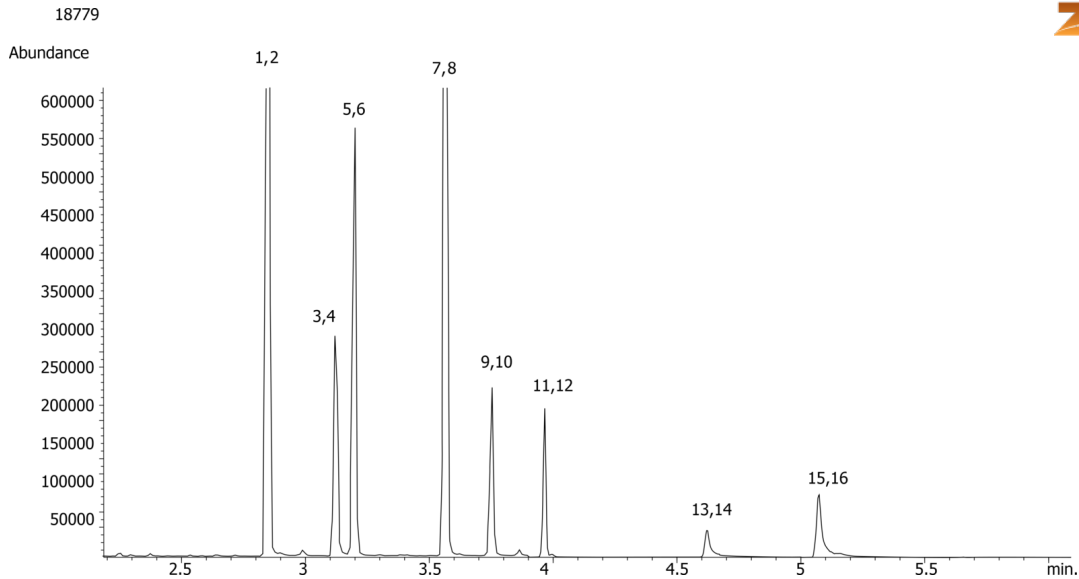


Benzodiazepines(with extra internal standards) by MSD on a ZB-Drug1**Column:** Zebtron™ ZB-Drug-1, GC Cap. Column 10 m x 0.18 mm x 0.18 µm, Ea**Phase:** Proprietary Drug Phase**Dimensions:** 10 meters x 0.18 mm x 0.18 µm**Order No:** 7CD-G023-08**Oven Profile:****Carrier Gas:** Constant Flow Helium, 0.7 mL/min**Injection:** Split 10.0000000000:1 0 µL @ 260°C**Detection:** Refractive Index (320°C)**Analyst Note:** Analytes are 4µg/mL derivatized with MTBSTFA

Sim ions:

(Quant ion*) D5-Nordiazepam (*332,334): Nordiazepam (*327,328,383): D5-Oxazepam (*462,519): Oxazepam (*457,513,459): D5-Temazepam (*362,288,390): Temazepam (*357,359,385): Lorazepam-d4 (*495,517)

**Products used in this application:**

Benzodiazepines(with extra internal standards) by MSD on a ZB-Drug1

ANALYTES:

- 1 Nordiazepam-d5
- 2 Nordiazepam
- 3 Diazepam-d5
- 4 Diazepam
- 5 Oxazepam-d5
- 6 Oxazepam
- 7 Lorazepam-d4
- 8 Lorazepam
- 9 Temazepam-d5
- 10 Temazepam
- 11 Clonazepam-d4
- 12 Clonazepam
- 13 Alprazolam-d5
- 14 Alprazolam
- 15 alpha-Hydroxyalprazolam-d5
- 16 alpha-Hydroxyalprazolam

