

EPA Method 552.2 - Haloacetic Acids on ZB-XLBHT

Column: Zebron™ ZB-XLB-HT, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea

Phase: Proprietary XLB Phase

Dimensions: 30 meters x 0.25 mm x 0.25 µm

Order No: 7HG-G024-11

Oven Profile:

Carrier Gas: Constant Flow Helium, 1.92 mL/min


Injection: Split 15.0000000000:1 0 µL @ 250°C

Detection: Refractive Index (310°C)

Analyst Note: Analytes are 250ppb in isooctane/MTBE



Products used in this application:

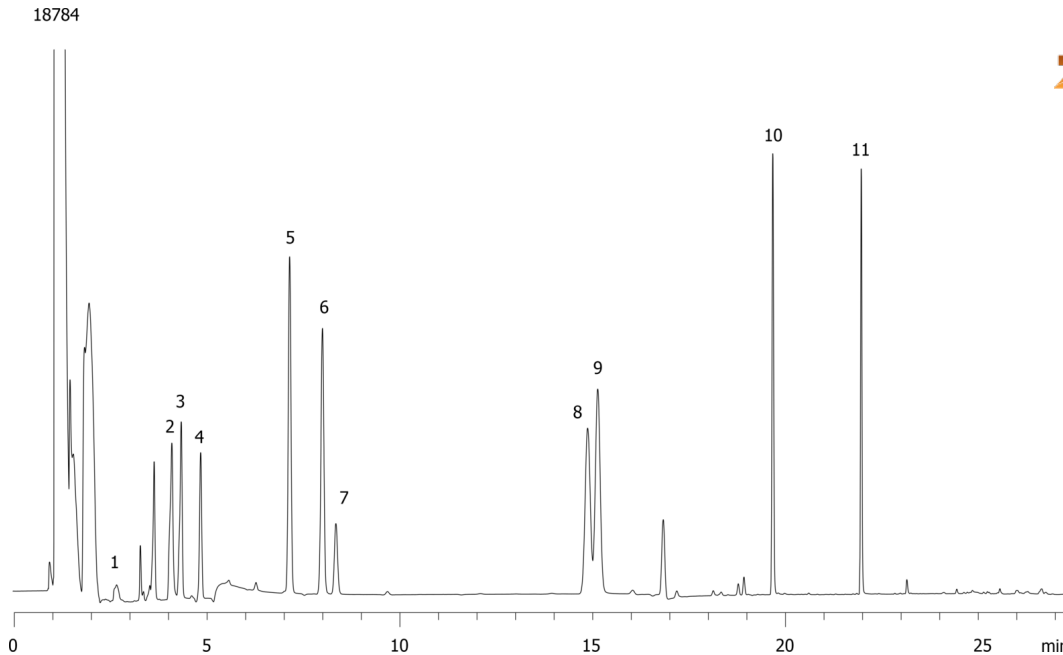

EPA Method 552.2 - Haloacetic Acids on ZB-XLBHT

ANALYTES:

- 1 Monochloroacetic acid
- 2 Monobromoacetic acid
- 3 Dichloroacetic acid
- 4 Dalapon
- 5 Trichloroacetic acid
- 6 Bromochloroacetic acid
- 7 1,2,3-Trichloropropane
- 8 Bromodichloroacetic acid
- 9 Dibromoacetic acid
- 10 Chlorodibromoacetic acid
- 11 Tribromoacetic acid

