


Atenolol and Related Substances CRS using Kinetex 2.6µm C18 100 x 4.6mm at 01.3ml/min**Column:** Kinetex® 2.6 µm C18 100 Å, LC Column 100 x 4.6 mm, Ea**Dimensions:** 100 x 4.6 mm ID**Order No:** 00D-4462-E0**Elution Type:** Isocratic**Eluent A:** EP Mobile Phase**Gradient** **Step No.** **Time (min)** **Pct A****Profile:** **1** 0 100**Flow Rate:** 1.3 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 226 nm (ambient)

18828

**Products used in this application:****ANALYTES:**

- 1** Impurity B
- 2** Impurity A
- 3** Impurity J
- 4** Impurity J
- 5** impurities D and E
- 6** Impurity F
- 7** Impurity G
- 8** Impurity H

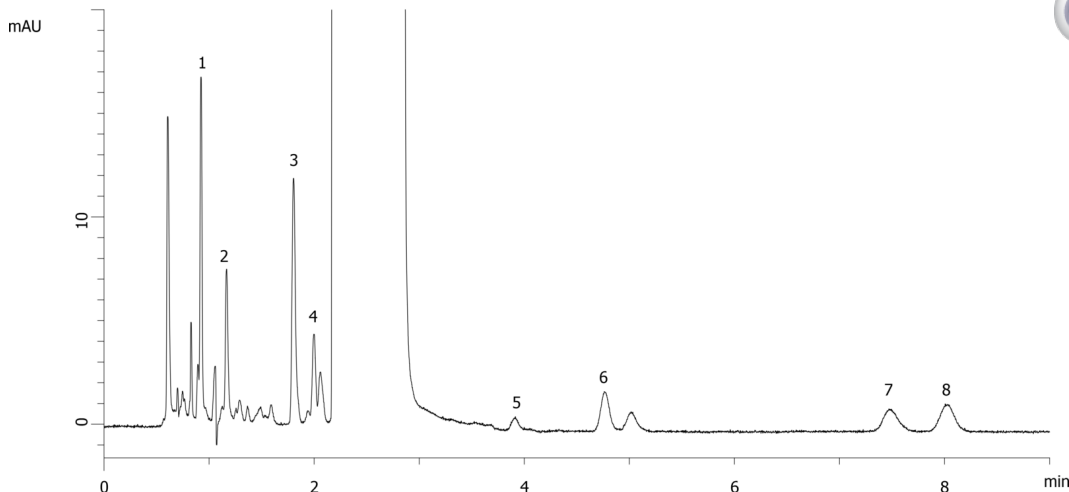