

Ornidazole on Lux 5 μ Cellulose-2 in NP 80% Hexane**Column:** Lux[®] 5 μ m Cellulose-2, LC Column 250 x 4.6 mm, Ea**Dimensions:** 250 x 4.6 mm ID**Order No:** 00G-4457-E0**Elution Type:** Isocratic**Eluent A:** 0.1% DEA in hexane**Eluent B:** 0.1% DEA in IPA

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	20	80	20

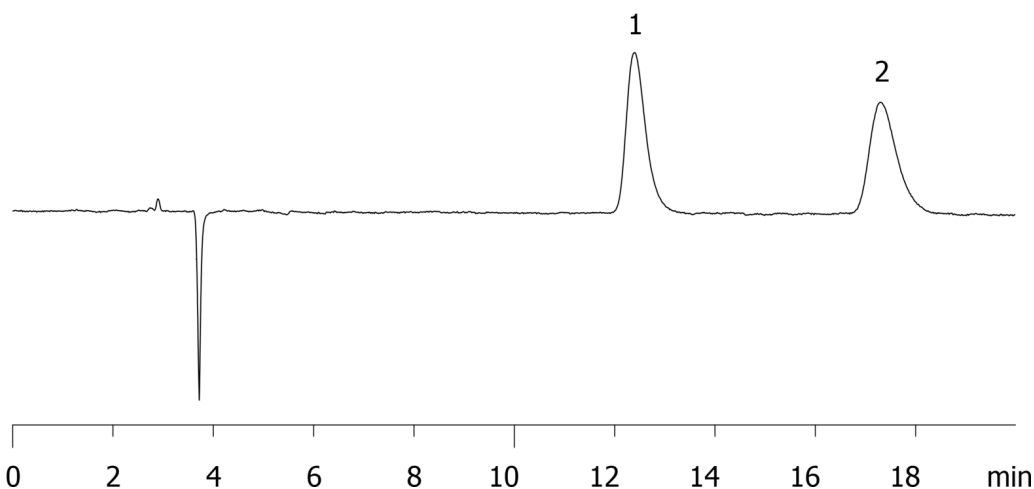
Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** Refractive Index @ 0.0000000000 (ambient)**Analyst Note:** SecurityGuard[™] Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Cellulose-2, 4 x 3.0mm ID, 10/Pk Part No.: AJ0-8366

- Holder Part No.: KJ0-4282

**Products used in this application:**

18831

**ANALYTES:****1** Ornidazole

Retention Time: 12.3 min

2 Ornidazole

Retention Time: 17.3 min

